

SPECIALIST LABORATORY SERVICES INFORMATION MANUAL

20th EDITION



	Page
TABLE OF CONTENTS	1
NATIONWIDE SPECIALIST LABORATORY	5
The Laboratory	5
Methodology	5
Quality Assurance	5
Personnel	6
Research and Clinical Trial Assays	6
Using our Services	6
Request Forms	6
Sample Collection Procedures	6
Sample Dispatch	6
Sample Results	6
Turnaround TimeSI units	7 7
Reference Ranges	7
Web Site	7
Specialist Clinical and Laboratory Endocrinology Advisors	8
About this manual	8
THYROID FUNCTION	9
Canine Hypothyroidism	10
Individual tests	10
Thyroid Panels and Dynamic Tests	12
Greyhound Thyroid Panels	13
Monitoring Therapy	14
Diagnosing Canine Hypothyroidism Flowchart	15
Feline Hyperthyroidism (Thyrotoxicosis)	16
Individual Tests	16
Thyroid Panels and Dynamic Tests	16
Feline Hypothyroidism	18
Equine Hypothyroidism	18

ADRENOCORTICAL FUNCTION	19
Canine Hyperadrenocorticism (HAC - Cushing's Disease)	20
Individual tests	20
Dynamic test protocols	20
Monitoring Therapy	24
 Monitoring Therapy Using "Pre-Pill" Single Cortisol Samples 	25
Diagnosing Canine Hyperadrenocorticism Flowchart	26
Canine Hypoadrenocorticism - Addison's Disease	27
Diagnosis	27
Monitoring therapy	29
Feline Hyperadrenocorticism (HAC - Cushing's Disease)	30
Guinea Pig Hyperadrenocorticism (HAC - Cushing's Disease)	32
Equine Pituitary Pars Intermedia Dysfunction (PPID)	33
Equine Metabolic Syndrome (EMS)	34
Adrenal Disease In Ferrets	35
Adrenal Sex Hormone Imbalance In Dogs (SHAP)	36
Canine Full Adrenal Sex Hormone Profile	37
Aldosterone	37
MALE REPRODUCTIVE FUNCTION	40
Cryptorchid/Testicular Function	41
FEMALE REPRODUCTIVE FUNCTION	44
Ovulation Detection	45
Ovarian Tissue Detection	46
Equine Inhibin	48
Equine AMH	48
PREGNANCY TESTS	50
Equine Progesterone (Early Indicative Pregnancy Test)	51
Equine Pregnant Mare Serum Gonadotrophin (PMSG)	51
Equine Oestrone Sulphate Pregnancy Test	51
Total Urinary Oestrogens (TUE)	51
Summary Table of Equine Pregnancy Tests	52
Canine Pregnancy Test - Relaxin	52
Feline Pregnancy Test - Relaxin	52
Faecal Pregnancy Detection	52
1 dood 1 regulatory Detection	52

GASTROINTESTINAL AND PANCREATIC	54
Trypsin-Like Immunoreactivity (TLI)	55
Lipase (DGGR)	55
Vitamin B12/Folate	56
Pancreatic Lipase Immunoreactivity (PLI)	56
Insulin	57
Fructosamine	58
CALCIUM DISORDERS	60
Parathyroid Hormone (PTH) and Ionised Calcium (iCa)	61
Parathyroid Hormone Related Peptide (PTHrP)	63
Vitamin D Metabolites	64
OTHER ENDOCRINE	66
Insulin-Like Growth Factor 1 (IGF-1)	67
Erythropoietin (EPO)	67
Gastrin (Special Request Assay)	67
Renin (Special Request Assay)	68
Glucagon (Special Request Assay)	70
Please note at the current time we are unable to perform Glucagon a withdrawal of the method we had validated. We will be looking at oth	•
NT-proBNP – Canine/Feline Cardiac Marker	71
Canine Benign Prostate Hyperplasia (BPH)	72
Normetanephrine/Creatinine Ratio for Phaeochromocytoma	72
THERAPEUTIC DRUG MONITORING	73
Phenobarbitone	74
Potassium Bromide	74
Levetiracetam (Keppra®)	74
Digoxin	75
Cyclosporine	75
Zonisamide	75

EN	DOCRINE THERAPEUTIC MONITORING	77
Bas	78	
Bas	78	
Vetoryl® Monitor Plus		78
'Pre Pill' Vetoryl® Monitor Plus		78
Basic Thyforon®		80
Thyforon® Monitor Plus		80
Basic Felimazole®		80
Fel	imazole® Monitor Plus	81
Fel	imazole® - diagnosis of feline hyperthyroidism	82
SE	ROLOGY/IMMUNOLOGY	83
Ace	etylcholine Receptor Antibodies (ACRAB)	84
Type 2M Muscle Fibre Antibody (2MAB)		84
Sarcoptic Mange Antibody (SMA)		84
Alle	ergy Testing	84
ON	ICOLOGY	85
Thy	ymidine Kinase	86
RE	FERENCES	87
AP	PENDICES	89
l.	Conversion Factors	90
II.	Assay Sample Requirements	91
III.	Specialist Biochemistry Profiles	92
IV.	Laboratory Request Form	93
V.	Laboratory Assay Schedule	94

Nationwide Specialist Laboratory

Nationwide Specialist Laboratories, has been established since 1998. Our lab was established by a small team of qualified scientists and technicians with many years of veterinary diagnostic and radioimmunoassay experience.

The aim of the laboratory is to find the best or "gold standard" assay technologies that can be made available to the veterinary profession at reasonable cost.

What we do

We provide unparalleled quality in all our endocrine assays giving the most accurate and precise results for every sample submitted thus ensuring you can confidently make the right diagnosis every time.

Our other laboratories

- NWL Poulton
- Abbey Veterinary Services

Full details can be found by visiting www.nwlabs.co.uk.

The Laboratory

Methodology

The measurement of hormones requires a higher level of analytical complexity than the "routine" clinical chemistries and usually requires techniques that depend on using antibodies to recognize the hormone molecules. These techniques are called immunoassays of which there are several kinds. Radioimmunoassay (RIA) techniques use low level radioactive materials as part of the detection system. This overcomes some of the problems of interference that can sometimes affect light-based detection systems such as enzyme immunoassay and chemiluminescence. Radioimmunoassay techniques are also more flexible and resilient when applied to a variety of species and sample types. Currently available chemiluminescence systems developed for human hormone analysis do not allow modification of assays to better suit veterinary patients. It is for these reasons that high quality radioimmunoassay techniques are the preferred choice of NationWide Specialist Laboratory.

Our interest in specialist methodologies has also given us the capabilities to expand our services to include special serology and drug monitoring.

Quality Assurance

We work to very high standards in laboratory practice to ensure you can have absolute confidence in our results. All analytical procedures and equipment are included in the internal quality system. All assay procedures are fully controlled using the relevant animal sera (Animal QC) and the laboratory participates in external quality assurance schemes where appropriate, including RIQAS for all biochemical analysis. All assays are fully validated for clinical use in every species if appropriate. NationWide Specialist Laboratory currently organizes and run the European Society of Veterinary Endocrinology (ESVE) External Quality Assessment Scheme.

Personnel

All NationWide Specialist Laboratory personnel are highly qualified technical officers with many years experience in veterinary endocrinology and analytical procedures. The diagnostic services of the laboratory are supported by access to several world-renowned veterinary clinical and laboratory endocrinologists.

Research and Clinical Trial Assays

We have a very wide range of knowledge, expertise and equipment enabling us to do virtually any type of esoteric assay, some of which may require complicated and sophisticated sample preparation or assay procedures. All of our staff are trained to very high standards in laboratory practice and can work and follow any required standards or protocols.

Using Our Services

Request Forms

Samples sent to the laboratory should be accompanied by an assay request form (obtainable from our Administration Department or from the web site (www.thehormonelab.com) alternatively details may be written on headed paper. To ensure the integrity of patient identification in the report, we request that the sample tube(s) be labelled very clearly and that dynamic tests have details of sample order or times as appropriate.

Sample Collection Procedures

Special sample preparation is not required for most of the tests but please refer to the request form, information sheets or the table at the back of this booklet for specific details. An overnight courier service is available for special samples (ACTH/PTH/PTHrP/Renin/Gastrin) for a small additional charge.

Sample Dispatch

Samples must be packed according to the current UN3373/P650 regulations for shipment of Diagnostic Specimens (a copy is available if required). Pre-paid Tracked Return labels for us can be printed from www.royalmail.com/track-my-return/create/4386 then just follow the on screen instructions.

Sample Results

All sample results are reported by fax and/or email and can be transferred direct to patient records. VetXML and FTP reporting options are also available. Urgent results can be obtained by telephoning the laboratory 8.00 am to 4.30 pm Monday to Friday. An answer machine operates out of hours for any messages.

Turnaround Time

Not all assays are run every day. Some are performed only twice weekly and certain low demand assays are run less often. The complexity of the high quality techniques we use can require long incubation times and numerous calibration and control samples necessitating the need to batch analyses to keep them affordable. An assay schedule is available on request but subject to change. Urgent requests can be analysed "off-schedule" but this will incur an additional charge (usually double list price).

SI units

It is the policy of the laboratory to report results in Systeme International (SI) units. A list of conversion factors for commonly measured analytes is provided as Appendix I.


Reference Ranges

A reference range list is available on request. Reference ranges are subject to revision if methods change; please refer to the sample test report for the current reference range(s).

Web Site

We have a very comprehensive web site which contains all of the information in this manual plus special information sheets and forms to download. It also contains news items and special promotions and offers so don't forget to keep looking!

www.thehormonelab.com

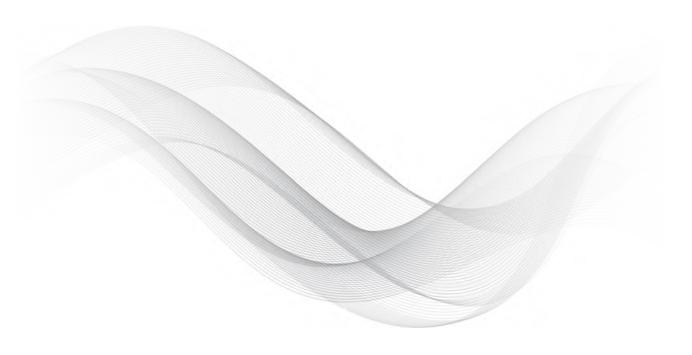

Specialist Clinical And Laboratory Endocrinology Advisors

Professor Michael E. Herrtage

M.A., B.V.Sc., D.V.Sc., D.V.R., D.V.D., D.S.A.M., F.R.C.V.S., Dip. E.C.V.I.M., Dip. E.C.V.D.I. RCVS and European Specialist in Small Animal Medicine



Mike Herrtage graduated from the Liverpool University and is now Emeritus Professor of Small Animal Medicine at the University of Cambridge and a Life Fellow of St. Edmund's College, Cambridge. He is a former Dean of the Cambridge Veterinary School and was in charge of the small animal medicine and diagnostic imaging services at the Queen's Veterinary School Hospital. His clinical interests include all aspects of small animal medicine and diagnostic imaging, but he has a particular interest in endocrine and metabolic disorders.

He was awarded the British Small Animal Veterinary Association (B.S.A.V.A.) Woodrow Award in 1986 for outstanding contributions in the field of small animal veterinary medicine and the B.S.A.V.A. Blaine Award for outstanding contributions to the advancement of small animal medicine in 2000. In 2014, he was awarded the World Small Animal Veterinary Association International Award for Scientific Achievement for outstanding contributions by a veterinarian, who has had a significant impact on the advancement of knowledge concerning the cause, detection, cure and/or control of disorders of companion animals. In 2016, he received the Kennel Club International Canine Awards' Lifetime Achievement Award in recognition of a career dedicated to improving canine health and welfare. He was voted Webinar Vet Speaker of the Year in 2015 and 2016 at the Annual Webinar Vet Veterinary Trust Awards. In 2019, he was awarded the B.S.A.V.A Bourgelat Award for really outstanding international contributions to the fields of small animal practice or science in their widest context by a qualified veterinarian.

He has been President of the British Veterinary Radiology Association, President of the British Small Animal Veterinary Association, President of the European Society of Veterinary Internal Medicine and President of the European Board of Veterinary Specialisation. He is a Diplomate of both the European College of Veterinary Internal Medicine and of the European College of Veterinary Diagnostic Imaging and was until recently President of the European College of Veterinary Internal Medicine. He was awarded the European Board of Veterinary Specialisation Award In 2018, for 'outstanding contributions to the course of Veterinary Specialisation in Europe'.

He has spoken at many international meetings and published over 200 articles in refereed journals.

Peter A Graham BVMS PhD CertVR

BVMS PhD CertVR DipECVCP MRCVS



Peter Graham graduated from

Glasgow University in 1989 where he remained as Small Animal House Physician and Research Scholar until 1996. During this period he was awarded the RCVS Certificate in Veterinary Radiology and a PhD on the Epidemiology and Management of Canine Diabetes Mellitus, Between 1996 and 2002 Peter was Assistant Professor at the world's largest specialist veterinary endocrinology laboratory in Michigan State University, USA, leading it as Section Chief from 2000. He was awarded Diplomate of the European College of Veterinary Clinical Pathologists in 2002. In 2002, he returned to the UK to take up the position of Managing Director of the NationWide Laboratories and NationWide Specialist Laboratories. In June 2014 he took up a position of Clinical Associate Professor Clinical Pathology and Endocrinology at the School of Veterinary Medicine and Science, The University of Nottingham.

He has a wealth of experience in the interpretation of endocrinology laboratory results, has authored and co-authored many publications and is a frequent invited speaker at international meetings.

About This Manual

This information manual describes the testing procedures available for a variety of endocrine diseases, pregnancy testing and therapeutic drug monitoring. The majority of the tests refer to canine, feline and equine submissions but we also have experience in other species. Please call if you have questions about the application of a measurement technique in a particular species.



Thyroid function



Thyroid Function

Canine Hypothyroidism

There are a great many factors affecting thyroid hormone levels in dogs. It is essential that a thorough clinical examination is performed and history taken prior to any diagnostic testing. The following factors are also very important to consider when carrying out thyroid function tests.

- Concurrent Therapy Certain therapies can affect thyroid hormone results and complicate their interpretation. Glucocorticoid and barbiturate medications often cause low total T4 (TT4) concentrations. Where possible, it helps to discontinue these therapies for 1 month prior to thyroid diagnostic testing. When clinical circumstances prevent withdrawal, a diagnostic panel including FT4ED (see below) should be selected. Sulphonamide products can cause a reversible hypothyroidism during their use and, consequently, thyroid diagnostic testing should be postponed until 3 weeks after they have been discontinued.
- <u>Breed Specific Ranges</u> Sight hounds are known to have lower TT4 levels compared to other breeds. The diagnosis of hypothyroidism in these breeds should be made cautiously and more reliance should be placed on other measures of thyroid function that are less influenced by breed such as TSH (or total T3).
- Low T4 State of Medical Illness (Sick Euthyroid Syndrome: SES) Dogs with non-thyroidal illnesses will often have low serum TT4 as part of their physiological response to that illness. This is not a T4 deficient state and thyroid supplementation is not appropriate. Instead, it is a reflection of the mechanism used to control metabolic rate during illness that is believed to improve the chances of survival. The effect of non-thyroidal illness on FT4ED values is less common and less dramatic. Ill dogs will often have depressed TT4 levels so the use of additional or alternative thyroid function tests along with careful evaluation of clinical signs and assessment of the likelihood of non-thyroidal illness are important in distinguishing the true hypothyroid dogs from those that are just responding to a non-thyroidal illness.

Individual tests

Total T4 (TT4)

Most dogs with hypothyroidism would be expected to have a low TT4 (~90% of them) making this a test with high, but slightly less than perfect, diagnostic sensitivity. On the other hand, the effects of non-thyroidal illness on TT4 means that many dogs with normal thyroid function will also have a low TT4 making the test poorly specific. Depending on the type of population sampled, up to 25% of dogs with normal thyroid function will yield a low TT4. This poor diagnostic specificity and less than perfect sensitivity means that TT4 has limited value as a stand-alone test for hypothyroidism. Diagnostic power is improved by combining it with TSH measurement or by performing a dynamic response test.

Total T4 also has a role to play in monitoring thyroid therapy where it can be used alone (if hours post-pill and dose frequency are recorded) or ideally in combination with TSH.

Free T4 by Equilibrium Dialysis (FT4ED)

Almost all (>99.9%) of circulating T4 is bound to carrier proteins leaving only a tiny fraction available to interact with tissues. This free fraction can be measured in an ultra sensitive radioimmunoassay following an equilibrium dialysis step. The analysis of FT4ED is the most accurate way of assessing the physiologically important thyroid status of an animal. Samples are dialysed, separating FT4 from serum proteins and protein bound T4. In most cases, TT4 and FT4ED will be highly correlated. The specific circumstances in which they are not are when we would recommend FT4ED as the thyroid hormone test of choice. It would be an advantage to measure FT4ED instead of, or in addition to, TT4 in the following situations:

- Non-thyroidal illness: one of the contributing mechanisms to the low TT4 we see in non-thyroidal illness is an alteration in thyroid hormone-protein binding. Although TT4 concentrations may be greatly reduced, the lower protein affinity for T4 means a higher fraction is available as free hormone and the FT4ED concentration usually remains within the reference range. This makes FT4ED a good test for distinguishing non-thyroidal illness from true hypothyroidism as the cause of a low TT4.
- Concurrent therapies: part of the effect on certain therapies on TT4 is mediated through thyroid hormone-protein binding meaning that FT4ED is less commonly and less dramatically affected by concurrent therapy. FT4ED is the analysis of choice when glucocorticoid or barbiturate therapies cannot be withdrawn prior to embarking on a thyroid diagnostic investigation.
- T4AA: the presence of T4 cross-reacting antibodies (T4AA) in the patient's serum will interfere with TT4 measurement causing false high values. The FT4ED procedure is unaffected by these antibodies because they are removed by the dialysis step. For an accurate estimation of thyroid status in a dog with T4AA, FT4ED is required. T4AA are present in the serum of approximately 10% of hypothyroid dogs as part of the thyroid pathologic process.

Canine TSH (cTSH)

We expect serum concentrations of cTSH to be high in animals with primary hypothyroidism because the negative feedback effect of thyroid hormones on pituitary production of TSH is lost. Indeed, this is the case most of the time but the diagnostic sensitivity is less than ideal. About 80-85% of hypothyroid dogs will have the expected high cTSH, unfortunately, leaving a proportion that will not. Conversely, cTSH measurement has good diagnostic specificity (up to 100%) meaning that false positives are rare. The combination of thyroid hormone analysis with cTSH measurement makes the most of the advantages of the individual tests while minimizing their deficiencies (see flow chart for interpretation of TT4 & cTSH).

Canine TSH may be measured at 30 minutes as part of a TRH stimulation test in the diagnosis of secondary hypothyroidism. Normal dogs should increase cTSH by at least 0.4 ng/mL.

Total T3 (TT3)

The analysis of TT3 is of little value in the diagnosis of hypothyroidism principally because of the high prevalence of cross-reacting T3 autoantibodies (T3AA) in hypothyroid dogs. These antibodies cause false results to be generated in T3 assays.

Thyroglobulin Autoantibody (TgAA)

The presence of TgAA in serum is strongly suggestive of immune mediated lymphocytic thyroiditis which is responsible for more than half of the cases of canine hypothyroidism. In the remainder there is no serological or histological evidence of inflammation. This test is used to document the presence and type of thyroid pathology. It does not provide information on the functional status of the thyroid glands. Evidence of thyroid function does not occur until lymphocytic thyroiditis has destroyed more than 50-60% of thyroid functional mass. Therefore TgAA evidence for thyroid pathology can be seen in animals before dysfunction occurs and while serum TT4 and cTSH concentrations are still normal. Some breeds of dog have a much higher prevalence of serum TgAA indicating a genetic predisposition to immune mediated thyroiditis. The screening of breeding lines or families for TgAA can be helpful for people wishing to breed away from a thyroiditis and hypothyroidism pre-disposition.

Anti T4/T3 Antibodies

Subsets of thyroglobulin antibody exist in certain hypothyroid dogs which cross-react with T4 or T3 assays. These antibodies, lead to falsely high thyroid hormone levels due to interference in the respective assay systems. FT4ED is not subject to this interference.

Thyroid Panels and Dynamic Tests

Basic Thyroid Profile (TP1)

Canine TSH is measured in combination with TT4 and this increases diagnostic accuracy over TT4 alone. This profile is the recommended single sample diagnostic screen.

Bronze Thyroid Profile (TP2)

Canine TSH is measured in addition to a TRH stimulation test. Dynamic function tests such as this have the potential to provide more information than single resting samples.

Copper Thyroid Profile (TP6)

TT4, Canine TSH and Free T4 ED. This panel provides the most comprehensive information available from a single sample at a more competitive price as TGAA is not measured and the possible cause of the thyroid pathology may not be considered important.

Silver Thyroid Profile (TP3)

Canine TSH is measured with Free T4 ED. This is a very useful combination for single sample analysis. This offers the advantages of TP4 but without the information on thyroid pathology provided by TgAA.

Gold Thyroid Profile (TP4)

Canine TSH and Free T4 ED are measured in addition to TgAA.

This panel does not suffer limitations of interference from non-thyroidal illness, concurrent therapy and T4AA that affect the measurement of TT4

Refer to Platinium profile (TP5) for suggested complete thyroid panel including TT4.

Platinium Thyroid Profile (TP5)

TT4, Canine TSH, Free T4 ED and TgAA. This panel provides the most comprehensive information available from a single sample. The TgAA identifies thyroid pathology and the TT4, FT4ED and TSH provide a complete picture of thyroid function.

TSH Stimulation Test

Unfortunately, pharmaceutical grade bovine TSH is no longer available. Chemical grade bovine TSH is available but extreme care should be taken if this product is used as there may be a risk of adverse reaction. A commercial pharmaceutical recombinant human TSH has been demonstrated to be useful for TSH stimulation testing in the dog. However, in its present form, it is very expensive.

- 1. Take blood sample for basal TT4 concentration.
- 2. Inject 0.1IU/kg TSH i/v. or 100 to 150 ug/dog human recombinant (rh)TSH
- 3. Take a second blood sample 4 6 hours later for post TT4 concentration.

Interpretation

T4 levels in a normal dog should increase by 1.5 - 2.0 times the basal concentration to reach a value above 26 nmol/L.

Greyhound Thyroid Panels

Greyhounds and other sighthounds (e.g. Saluki's) have total serum thyroxine (TT4) levels that are lower than those of other breeds of dog. Thyroid investigation is not uncommon in greyhounds as part of the work up for poor performance, behavioural change, bald thighs and fertility concerns. NationWide Specialist Laboratory have put together profiles of thyroid tests specifically for the greyhound breed.

Greyhound Thyroid Investigation (GTH1)

Sensitive TT4, TT3, TSH and TGAA

Greyhound Thyroid Monitoring (GTH2)

Sensitive TT4, TT3 and Canine TSH.

TRH Stimulation Test

- 1. Take blood sample for basal TT4 concentration.
- 2. Inject TRH i/v slowly over one minute.
- 3. 1 5 kg 100ug TRH 5 - 30 kg 200ug TRH >30 kg 300ug TRH
- 4. Take a second blood sample 4 6 hours later for post TT4 concentration.

Interpretation

TT4 levels in a normal dog should increase by about 1.2 times the basal concentration to reach a value above 25 nmol/L. If the pre stimulation sample is above 25 nmol/L the dog is

Thyroid Function

likely to be normal, regardless of the post stimulation TT4.

Hypothyroid dogs usually show low basal TT4 levels, which fail to respond to TRH or stimulate to a value below 25 nmol/L. Sometimes a high/normal basal concentration may be seen which fails to increase 1.2 times. On these occasions it is likely that the thyroid gland is being stimulated maximally and cannot respond further (in the absence of T4AA).

Diagnosis of secondary hypothyroidism can be attempted by measuring cTSH at zero and 30 minutes post TRH stimulation. Normal dogs should increase cTSH by at least 0.4 ng/mL.


Monitoring Therapy

Thyforon® (Dechra Veterinary Products), Thyroxanil® (Animalcare) and Leventa® (MSD) are the only licensed veterinary preparations of levothyroxine. The initial recommended dose varies by datasheet but is around 20ug/kg daily. However, a large proportion of hypothyroid dogs can be successfully treated with a lower dose. TT4 levels should be monitored after several weeks of therapy to determine whether a dose change would be appropriate. For animals receiving therapy twice daily, the time post pill at which the monitoring sample is taken is less important than that it is recorded.

TT4 levels should ideally be between 30 - 47 nmol/L at around the time of expected peak concentrations (3 hours post pill) and be above 19.0 nmol/L by the time the next pill is due (trough concentration). The half-life of TT4 in the dog can vary but is in the region of 10 to 14 hours.

The measurement of both TT4 and cTSH is recommended in monitoring thyroid therapy in dogs. cTSH reflects the adequacy of thyroid replacement therapy in the preceding days, not just the day of the test. This can help identify inconsistencies in dosing and also prevent inappropriate management decisions being made based on unrepresentative single day TT4 results.

Diagnosing Canine Hypothyroidism Flowchart

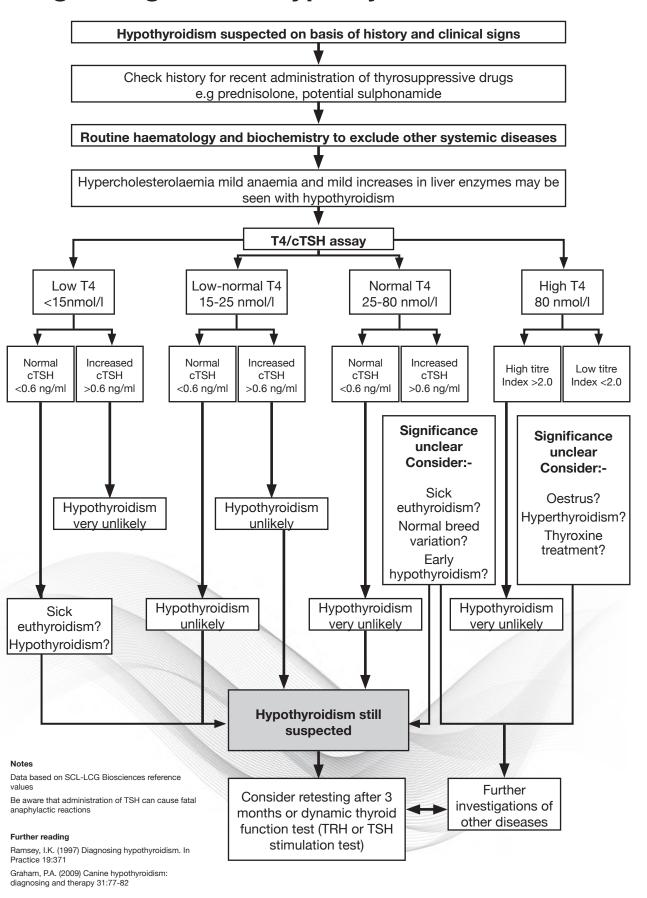

Feline Hyperthyroidism (Thyrotoxicosis)

Individual Tests

Total T4

TT4 is generally used as the main diagnostic test for hyperthyroidism and to monitor T4 levels post treatment. TT4 should be used as the first screening test. If TT4 is raised (> 60 nmol/L) it is very likely the cat is hyperthyroid. If TT4 is within the reference range then further tests may be helpful in ruling out or confirming hyperthyroidism. Often hyperthyroid cats with concurrent non-thyroidal illness will have normal TT4 concentrations (occult hyperthyroidism). In these cases the detection of a palpable thyroid nodule will be a strong indicator that hyperthyroidism should still be suspected and investigated with FT4ED or a dynamic function test. If TT4 is in the lower end of the reference range (<30 nmol/L) it is unlikely that the cat is hyperthyroid. Euthyroid cats with significant non-thyroidal illness will have low-normal or subnormal TT4 concentrations; generally, the more severe the illness the lower the TT4.

Free T4 by Equilibrium Dialysis (FT4ED)

FT4ED is the most accurate way to measure free T4 and has greater diagnostic sensitivity for hyperthyroidism in sick animals where the TT4 may be depressed into the reference range. FT4ED is less affected by altered binding protein characteristics seen in non-thyroidal illnesses as the samples are dialysed prior to assay. FT4ED levels are elevated in hyperthyroid cats that have TT4 levels within the normal range (19 - 65 nmol/L). Care is needed when using FT4ED as a stand-alone test because it may be elevated in sick euthyroid cats. To mitigate this risk of false positives, FT4ED should be measured in conjunction with TT4 and also when there are appropriate clinical signs such as a palpable thyroid nodule.

Thyrotropin (TSH)

Using the same assay as is used to measure canine TSH there is sufficient cross-reactivity for us to detect increased levels of feline TSH. This can be a helpful analyte in the identification of iatrogenic hypothyroidism following treatment of hyperthyroidism and in the less common presentation of naturally occurring feline hypothyroidism (both congenital +/- goitre and inflammatory forms have been recognised). The combination of azotaemia and iatrogenic hypothyroidism has been shown have a detrimental effect on the survival of treated hyperthyroid cats.

Thyroid Panels and Dynamic Tests

Feline Hyperthyroid Profile

The combination of TT4 and FT4ED offers the advantages of improved diagnostic sensitivity (less false negatives) over TT4 alone and improved diagnostic specificity (less false positives) compared to FT4ED alone.

Feline Hyperthyroid Profile - Silver

(Total T4, ALT, ALP, Total Protein, Urea, Creatinine, Phosphorus). Specific biochemistry tests in addition to TT4.

Feline Hyperthyroid Profile - Gold

(As silver plus FT4ED to improve diagnostic sensitivity and specificity).

T3 Suppression Test

This is a useful test for diagnosing borderline cases of hyperthyroidism where the TT4 is consistently within the normal range. It is very important to measure Total T3 (TT3) as well to ensure the cat actually absorbed the T3 dose.

- 1. Collect basal blood sample for TT4 and TT3. Store separated serum frozen until protocol completed and submit with post-sample
- 2. Administer T3 (Tertroxin) orally every 8 hours for a total of 7 doses according to body weight:
 - cats <5 kg 20 micrograms of Tertroxin cats >5 kg 30 micrograms of Tertroxin
- 3. Collect second blood sample 2 6 hours after final dose
- 4. Label both samples clearly and request TT4 and TT3.

Interpretation

TT3 suppression normal cats usually show at least 50% reduction in TT4 levels following suppression. Hyperthyroid cats generally show limited suppression. TT3 should be increase indicating the T3 tablets have been absorbed.

TRH Stimulation Test

Can be useful in cases of mild or borderline hyperthyroidism.

- 1. Collect basal blood sample.
- 2. Administer 100 µg/kg i/v slowly over one minute
- 3. Collect second blood sample 4 hours later
- 4. Label samples clearly and request TT4.

Interpretation

Normal cats increase 1.5 - 2 times following stimulation.

Hyperthyroid cats generally show little or no rise in TT4 levels from a high or high-normal baseline.



Feline Hypothyroidism

Feline hypothyroidism is considered to be a rare condition. Options for diagnosis include the measurement of serum TT4 or Free T4 and TSH (TP1 or TP3 – see canine hypothyroidism), or a TRH stimulation test. A TSH stimulation test would be ideal if a TSH preparation is available (including recombinant human TSH; Thyrogen®, Genzyme; 0.025 to 0.200 mg IV of rhTSH).

The most common form of feline hypothyroidism is iatrogenic following treatment of hyperthyroidism. Less commonly we recognise naturally occurring feline hypothyroidism (both congenital +/- goitre and inflammatory forms exist). The combination of azotaemia and iatrogenic hypothyroidism has been shown to have a detrimental effect on the survival of treated hyperthyroid cats.

TRH Stimulation Test

- 1. Collect blood for basal TT4 concentration
- 2. Inject 100 µg TRH i/v slowly over one minute. Side effects including salivation and vomiting can be seen occasionally
- 3. Collect second blood sample 4 6 hours later
- 4. Label samples clearly and request TT4

Interpretation

Normal cats - post stimulation result is expected to be about 1.5 - 2 times greater than the basal T4 concentration.

Monitoring of feline hyperthyroidism is detailed at the end of the booklet under endocrine therapeutic monitoring.

Equine Hypothyroidism

The prevalence of equine hypothyroidism is controversial. Most believe it to be a rare condition.

TRH Stimulation Test

- 1. Take a basal blood sample
- 2. Inject 1 mg TRH per horse or 0.5 mg per pony i/v slowly over one minute
- 3. Collect a second blood sample 4 6 hours later
- 4. Label both samples clearly and request TT4 analysis.

Interpretation

In normal horses and ponies, TT4 is expected to rise at least 1.4 times basal following stimulation.

In the horse, TT3 rises earlier and more dramatically than TT4 and so its measurement may be valuable. In normal animals TT3 will rise up to 3 times the basal value after 4 - 6 hours.



Adrenocortical Function



Adrenocortical Function

Canine Hyperadrenocorticism (HAC - Cushing's Disease)

Hyperadrenocorticism occurs in middle–aged and geriatric dogs. In around 80% of cases it is pituitary-dependent (PDH). In the remainder, the disease results from primary adrenal-dependent neoplasia or hyperplasia (ADH).

Individual tests

Cortisol

Basal cortisol levels are of no value in diagnosing HAC. There is significant overlap in cortisol concentration among healthy, stressed, sick and affected dogs.

Cortisol/Creatinine Ratio (CCR)

This is a very sensitive test to exclude HAC (very few false negatives) but must not be used to diagnose HAC as it is very poorly specific and non-adrenal illness commonly gives a positive result. A morning urine sample is collected in the animal's home environment. This reflects the cortisol released over several hours. Normal dogs will have a CCR less than 30×10^{-6} .

Dynamic test protocols

Which test should be used to diagnose HAC reliably in the dog?

Unfortunately, there is no perfect diagnostic test or test protocol for hyperadrenocorticism. The line between physiologic adrenal stress responses to non-adrenal illness and pathologic hyperfunction of the pituitary-adrenal axis can be fine indeed. As a consequence, false positive and false negative test results occur. The different available diagnostic protocols available have different properties, advantages and disadvantages but none is perfect. In establishing confidence that hyperadrenocorticism is confirmed or ruled out more than one diagnostic test protocol may be required.

ACTH stimulation is a good screening test in the first instance and is the test of choice for diagnosing iatrogenic Cushing's and monitoring anti-adrenal therapy. It has a lower false positive rate than the low-dose dexamethasone suppression test but a significant false negative rate. It will reliably diagnose about 85% of PDH cases but only 50% of ADH cases. It is quick and simple to perform and is less affected by stress and non-adrenal illness. The initial values are useful as a reference to monitor effectiveness of treatment.

Low-dose Dexamethasone screening test is more sensitive than the ACTH stimulation test in confirming HAC. Diagnostic sensitivity is almost perfect (>95%) meaning that false negative results very seldom occur giving us tremendous confidence in a negative test result. Unfortunately, false positive results are common especially when there is concurrent non-adrenal illness or other sources of stress. Positive results of the low-dose dexamethasone test should be regarded with suspicion in dogs known to have significant non-adrenal illness. Ideally, the test should be postponed until any identified non-adrenal illnesses have been resolved or stabilised. See "Diagnosing Canine Hyperadrenocorticism" flowchart (page 27).

17-hydroxyprogesterone measurement in the diagnosis of adrenocortical disease


Certain functional adrenocortical tumours do not have cortisol as their principle secretory product but instead the tumour has developed in such a way that it is adrenal steroid precursors that are released into the circulation. Many dogs with this type of functional adrenal tumour will have clinical signs suggestive of hyperadrenocorticism. A common cortisol response to ACTH in these cases is a "flat-line, mid-range" pattern. Measuring 17OHP before and after ACTH stimulation (same sample as you would take for cortisol) can be very helpful in identifying these tumours.

Additionally, in cases suspected of hyperadrenocorticism which do demonstrate an exaggerated cortisol response to ACTH, the additional measurement of 17OHP can help improve the confidence in ruling hyperadrenocorticism in or out. Less than 10% of dogs with classic hyperadrenocorticism would be expected to have post-ACTH 17OHP <4.5 nmol/L (Chapman et al, 2003 Veterinary Record, 153(25):771-5) and fewer than 10% of dogs which do not have classic hyperadrenocorticism would be expected to have post-ACTH 17OHP >16.7 nmol/L.

Consider measuring 170HP in addition to cortisol when:

- Flat-line, mid-range cortisol response to ACTH
 - Investigate functional adrenocortical tumour
- Mildly exaggerated cortisol response to ACTH and few appropriate clinical signs of HAC
 - Look for post-ACTH <4.5nmol/L to improve rule-out
- Normal cortisol response to ACTH and clinical signs consistent with HAC
 - Investigate "atypical" hyperadrenocorticism (Ristic et al 2002 JVIM 16(4):433-9) and "Alopecia X" (page 35)

170HP may also be useful in ruling out HAC where an exaggerated Cortisol response is seen which is mainly due to stress, often associated with other illness. The Cortisol levels may be borderline or high but the 170HP response to ACTH is normal.

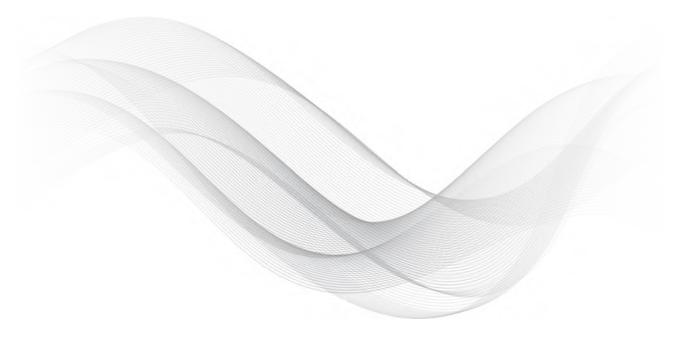

ACTH Stimulation Test

- 1. Collect basal blood sample.
- 2. Immediately inject 0.25 mg or 5ug/kg synthetic ACTH (Synacthen /Tetracosactide) i/v or i/m
- 3. Collect a second blood sample one hour later.
- 4. Label samples clearly and request Cortisol analysis.

Interpretation

Normal dogs will show an increase in cortisol levels of up to 450 nmol/L post stimulation.

An exaggerated response is expected in animals with PDH and cortisol concentrations rise above 600 nmol/L and often above 1000 nmol/L.

Dogs with an adrenal tumour often have basal cortisol above 250 nmol/L with little or no change after stimulation. However the ACTH stimulation test is not as sensitive at detecting adrenal tumours and negative results should be confirmed with a low dose dexamethasone test when strong clinical suspicion remains.

Low Dose Dexamethasone Screening Test

- 1. Collect a basal blood sample
- 2. Inject 0.01 mg/kg i/v of Dexamethasone.
- 3. Collect two further blood samples 3 hours and 8 hours later
- 4. Label sample times clearly on the tubes and request cortisol.

Interpretation

The low-dose dexamethasone screening test is interpreted in two stages. Firstly, the presence or absence of hyperadrenocorticism is established by examining the 8-hour result. A value above 40 nmol/L is positive. The next step applies in the positive cases and checks for evidence of suppression sufficient to identify the source of hyperadrenocorticism. Greater than 50% suppression from the baseline value at either 3 or 8 hours is consistent with pituitary dependent disease. When there is no, or minimal, evidence of suppression it is necessary to follow-up with a differentiation test to identify the source (high-dose dexamethasone suppression test, endogenous (plasma) ACTH, adrenal ultrasonography). Both dogs with an adrenal tumour and some dexamethasone "resistant" pituitary cases will not depress by 50%. In about 60% of hyperadrenocorticism cases, we will obtain differentiation from the low-dose test alone with no need to move on to the differentiation tests.

High Dose Dexamethasone Suppression Test

This test is used to help distinguish between PDH and an adrenal tumour. It must only be used after Cushing's disease has been confirmed by other tests (low dose dex. or ACTH stim.)

The disadvantages to the HDDST are the 8 hour duration and that a small proportion of pituitary dependent cases will fail to suppress even at substantial doses of dexamethasone. An alternative differentiation test is the endogenous (plasma) ACTH which has the advantage of being a single point test.

- 1. Collect a basal blood sample.
- 2. Inject 0.1 1.0 mg/kg i/v of Dexamethasone.
- 3. Collect two further blood samples 3 hours and 8 hours later
- 4. Label sample times clearly on the tubes and request cortisol.

Interpretation

Any suppression greater than 50% of the baseline concentration indicates a pituitary source. A failure to suppress by 50% is consistent with either an adrenal tumour or a "dexamethasone-resistant" pituitary lesion. When no suppression is observed, repeating the test with a higher dose of dexamethasone or diagnostic imaging of the adrenals and/or pituitary could be considered.

Endogenous (plasma) ACTH

Endogenous ACTH may be used as a reliable tool for differentiating PDH from ADH when hyperadrenocorticism has already been confirmed by ACTH stimulation test or low-dose dexamethasone screening test. However, it has no value as a diagnostic test for hyperadrenocorticism in dogs and cats. Special sample handling procedures apply.

ACTH Sample Preparation

Sample MUST be EDTA plasma or EDTA plasma with Aprotinin.

- 1. Request a transport pack and follow the instructions that come with the pack for sample preparation.
- 2. Take the blood sample into the **cooled Aprotinin EDTA tube** (supplied in pack) which has been cooled in the fridge.
- 3. Mix very well but gently at least 10 times and centrifuge as quickly as possible.
- 4. Transfer the plasma into a cooled plastic (not glass!) PLAIN tube kept in the fridge.
- 5. Immediately freeze (<-10°C) the plasma sample and keep frozen until dispatch in the transport pack. Freeze the sample first, **separately**, before placing in the frozen pack for transport.
- 6. Mark the tube with "EDTA plasma" and animal/owner name.

Cortisol/Creatinine Ratio (CCR) with Dexamethasone Suppression Test

When Cushing's disease is strongly suspected, this test can be used when the dog is easily stressed or difficult to sample. The owner does all the sampling at home.

- 1. Day 1, collect a first morning urine sample.
- 2. Mix the urine and add some to Sample Tube 1, place in the fridge until dispatch.
- 3. Day 2, collect a first morning urine sample.
- 4. Mix the urine and add some to Sample Tube 2, place in the fridge until dispatch.
- 5. Immediately after urine collection, note the time and give the dog the required number of Dexamethasone tablets (dose = 0.1mg/kg).
- 6. 8 hours later give the dog a 2nd set of Dexamethasone tablets.
- 7. 16 hours later give the dog a 3rd set of Dexamethasone tablets.
- 8. Day 3, collect a first morning urine sample.
- 9. Mix the urine and add some to Sample Tube 3, place in the fridge until dispatch.
- 10. Send all 3 urine samples to the laboratory for CCR.

Interpretation

HAC is suspected if the CCR is greater than 30×10^{-6} in two consecutive morning urine samples.

If the CCR in the 3rd urine sample is depressed by 50% of the mean CCR of the first two samples, PDH is the likely diagnosis.

If the suppression is less than 50%, an endogenous ACTH is suggested to confirm an adrenal tumour

Normal dogs will have a CCR less than 30 x 10⁻⁶.


Monitoring Therapy

The standard ACTH stimulation test is the test of choice in monitoring anti-adrenal therapy (Lysodren, Vetoryl'). The aim of therapy is to achieve a modestly sub-normal ACTH response in order to abolish clinical signs.


Dogs receiving Vetoryl must have the ACTH response test done 4 – 6 hours after the Vetoryl tablet. The post ACTH cortisol level should be between 40 and 200 nmol/L to reflect good control. Experience in the field has shown that a post ACTH cortisol of up to 250 nmol/L may not require an increase in dosage for continued clinical effect. It is important that there is some stimulation from the pre to the post sample. Cortisol levels of <20 (pre) and <20 (post) are not acceptable. Clinical signs must be monitored closely. ACTH stimulation tests should be carried out at 10 days, 4 weeks, 12 weeks and then every 3 months post treatment with Vetoryl.

Monitoring Therapy Using "Pre-Pill" Single Cortisol Samples

The standard ACTH stimulation test is the documented procedure in the data sheet to monitor Vetoryl treatment, however research carried out by Professor Ian Ramsey at Glasgow Vet School has shown that "pre-pill" single cortisol samples may be an effective alternative to ACTH stimulation tests when monitoring Vetoryl when used in conjunction with the animals condition (well or unwell) and any clinical signs of HAC. Advice regarding pre-pill single cortisol testing may be obtained from <a href="mailto:test-needle-bracked-less-signs-ne

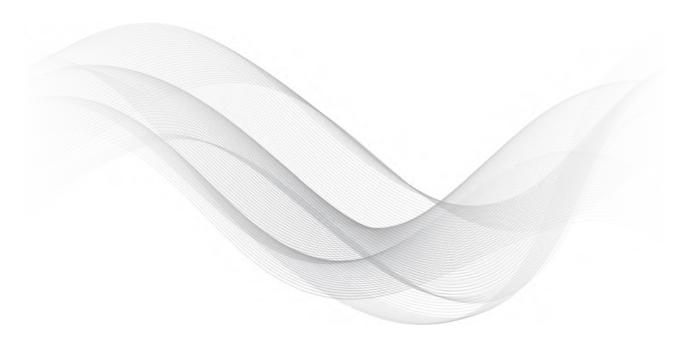

Diagnosing Canine Hyperadrenocorticism Flowchart

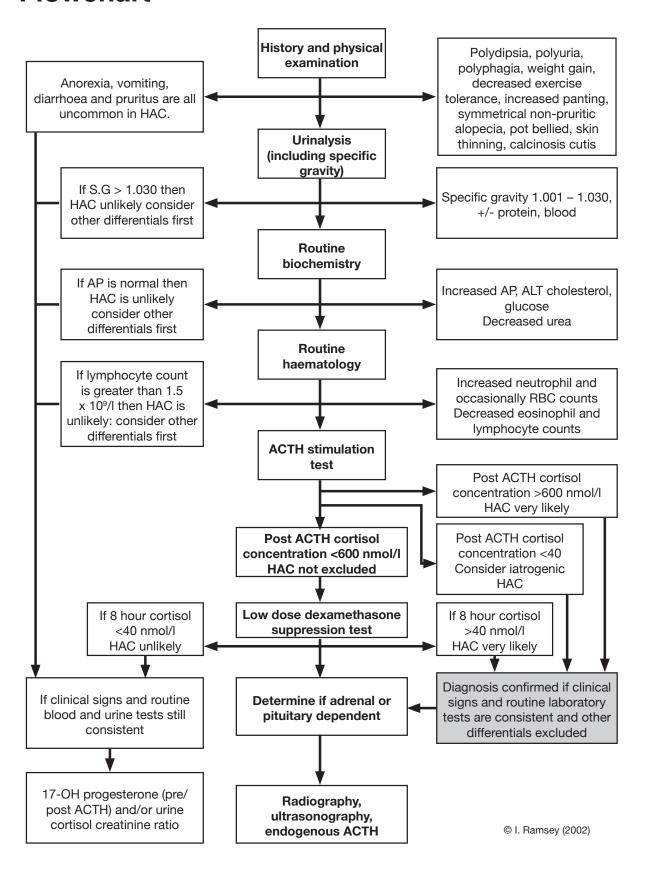

Canine Hypoadrenocorticism - Addison's Disease

Diagnosis

Primary hypoadrenocorticism results when more than 90% of both adrenal cortices are destroyed and this leads to a clinical deficiency of all adrenal hormones. The ACTH stimulation test will provide confirmation of Addison's. It is also useful to measure the Sodium/Potassium ratio as evidence for mineralocorticoid deficiency. A low Na:K ratio may not be seen in some cases of primary hypoadrenocorticism until a very advanced stage of the disease. That is, an animal can have confirmed primary hypoadrenocorticism without abnormal electrolytes. A diagnosis of "atypical" or "secondary" hypoadrenocorticism can only be made if failure to produce aldosterone (see below) is confirmed in addition to failure to produce cortisol.

ACTH Stimulation Test

- 1. Collect basal blood sample.
- 2. Immediately inject 0.25 mg or 5ug/kg synthetic ACTH (Synacthen/Tetracosactide) i/v or i/m
- 3. Collect a second blood sample one hour later.
- 4. Label samples clearly and request Cortisol analysis.

Interpretation

The basal cortisol will be low with little or no response after ACTH. A classic Addison's ACTH response test would be Pre cortisol = <20 nmol/L, Post cortisol = <20 nmol/L. False positive diagnoses occur in animals which have been receiving glucocorticoid therapy recently including oral, injectable and topical (eye, ear and skin) preparations. For confirmation of primary hypoadrenocorticism in dogs which have been receiving glucocorticoid therapy, consider measuring aldosterone instead of, or in addition to, cortisol.

Prednisolone, prednisone and hydrocortisone will cross-react in the cortisol assay and give falsely elevated results. Such therapies should be avoided in the 24 to 36 hours prior to the test. Dexamethasone does not have this effect so it can safely be used prior to confirmation of diagnosis without causing assay interference. However, prolonged use (several days) of dexamethasone prior to performing the ACTH response test can increase the risk of false positive diagnosis and require the measurement of aldosterone for confirmation.

Endogenous (plasma) ACTH

ACTH concentrations can be useful in distinguishing primary from secondary hypoadrenocorticism (see differentiation in Canine hyperadrenocorticism section). Dogs with primary Addison's have very high levels of ACTH, usually greater than 500 pg/mL. Dogs with secondary Addison's have low or undetectable (< 5.0 pg/mL) levels of ACTH.

Aldosterone

Aldosterone measurement may also be used to differentiate between primary and secondary Addison's and this can be performed on the same samples used to measure the cortisol, with no difference in sample handling procedures.

In primary hypoadrenocorticism, there will be failure of both glucocorticoid and mineralocorticoid production. In secondary hypoadrenocorticism (pituitary ACTH production failure, exogenous steroids) glucocorticoid production will have failed but mineralocorticoid production is unaffected. The trophic support for mineralocorticoid production comes from the renin-angiotensin system not ACTH. However, there is a sufficient aldosterone response to ACTH stimulation to use this as a dynamic diagnostic test for mineralocorticoid production which has advantages over measuring baseline aldosterone alone.

Measuring aldosterone in addition to cortisol is important when it is necessary to perform an ACTH stimulation test to confirm a diagnosis of hypoadrenocorticism after an animal has already received symptomatic glucocorticoid therapy. It is also essential if a diagnosis of atypical or secondary hypoadrenocorticism is being considered (unless cACTH is measured).

Cortisol-to-ACTH Ratio (CAR) for Canine Hypoadrenocortism (HA).

When there is a problem with the supply of tetracosactide (Synacthen). The Cortisol to ACTH ratio (CAR) can provide clinicians with a reliable alternative test to diagnose hypoadrenocortism (HA) (Addison's disease).

A publication in the JVIM looked at the CAR in dogs with HA, healthy dogs and those with diseases mimicking HA. (Boretti and others 2015).

Measurement of baseline Cortisol has been shown to be a useful screening test to rule out HA as the disease is very unlikely with a baseline Cortisol > 55 nmol/L but the low specificity of this test means an ACTH stimulation test must be performed to safely exclude HA (Bovens and others 2014, Lennon and others 2007).

Plasma ACTH concentrations are significantly higher in dogs with HA compared to dogs with diseases mimicking HA and in healthy dogs. The CAR is significantly lower in dogs with HA compared to dogs with diseases mimicking HA and in healthy dogs (Boretti and others 2015, Javadi and others 2006).

A diagnostic sensitivity of 100% and a specificity of 99% has been reported with a CAR of > 0.003 nmol/pg (when measured using cortisol nmol/l and ACTH in pg/ml). Values less than 0.003 nmol/pg strongly support a diagnosis of HA – most dogs with HA are less than 0.0003 nmol/pg (Boretti and others 2015). Values above 0.003 suggest it is highly unlikely that the patient has HA, most dogs with diseases that mimic HA are above 0.08 nmol/pg (Boretti and others 2015). This test also works with so- called 'atypical' hypoadrenocorticism (Lathan and others 2014). It should be noted however that the exact cut-off has varied between the 3 studies so far published and caution is required near to the cut-off value (Boretti and others 2015). Current steroid administration must be excluded.

The careful collection and handling of samples for the measurement of plasma ACTH sample is **critical** and special procedures must be followed to avoid any loss of ACTH activity or this may make the results misleading. Briefly, a single blood sample is collected into EDTA and then the plasma is rapidly separated using a centrifuge. The plasma is then taken into a serum tube and placed in a freezer. From collection to freezing must be less 15 minutes or less. The sample must be sent by a guaranteed express delivery service in a special freezer pack obtained from NationWide Specialist Laboratories.

Adrenocortical Function

Please telephone for further advice on sampling and ordering the special freezer pack.

NationWide Specialist Laboratories have a very sensitive canine ACTH assay and are able to offer the analysis and interpretation of CAR to help in the diagnosis of HA.

Canine ACTH is normally run on Monday but we are prepared to run same day assays if required, please ask for details.


Monitoring therapy

The laboratory monitoring of primary hypoadrenocorticism focuses on electrolyte status. There is no value in repeating ACTH stimulation tests in animals which have failed adrenal glands. The artefactual cortisol concentrations generated from the cross-reaction of prednisolone in the assay are not biologically equivalent and so a cortisol assay cannot assess the adequacy of glucocorticoid therapy. Biologically equivalent cortisol results can be obtained only during periods that hydrocortisone is the exclusive corticosteroid therapy.

The monitoring of secondary hypoadrenocorticism is based mostly on clinical response with periodic attempts to find the lowest glucocorticoid dose at which the animal appears well. Liver enzyme analysis can help identify steroid hepatopathy resulting from glucocorticoid excess.

Zycortal® (Dechra Veterinary Products) is indicated for replacement therapy for mineralocorticoid deficiency in dogs with primary hypoadrenocorticism. Zycortal contains desoxycortone pivalate (DOCP) – an aldosterone analogue with purely mineralocorticoid effect. Zycortal is available as a prolonged-release suspension for subcutaneous injection and the dosing interval is approximately 25 days. The exact dose to be given and the interval between treatments depend on the dog's clinical response as well as sodium and potassium levels. A glucocorticoid (e.g. prednisolone) is also needed daily to replace the missing cortisol.

To obtain a copy of the Zycortal[®] monitoring and dose adjustment guidelines (5 Golden Rules), or for further advice, contct Dechra by emailing **technical@dechra.com**.

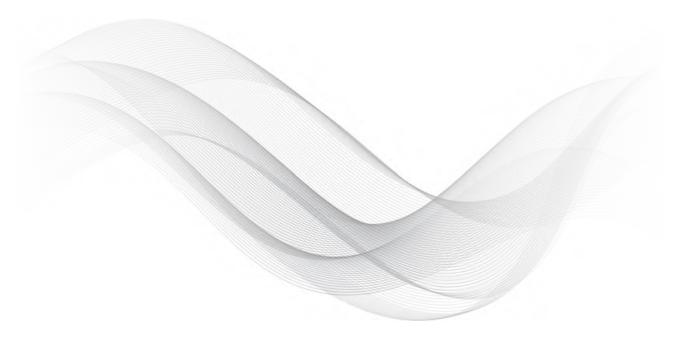

Feline Hyperadrenocorticism (HAC - Cushing's Disease)

Hyperadrenocorticism is a rare condition in the cat, much more so than in the dog. The scarcity of true clinical cases means that there are few studies of diagnostic accuracy for the different screening tests. Unlike the situation in the dog, we do not have good information on diagnostic sensitivity and specificity. In its absence, recommendations have been derived from studies in healthy cats and some extrapolated from the dog.

The clinical appearance of reported cases has been dramatic. Alopecia has not been a feature but potbelly and dramatic skin thinning have been. The skin may become so thin that it tears spontaneously or on handling.

The cat does not have a steroid induced iso-enzyme of alkaline phosphatase like that in the dog.

If feline hyperadrenocorticism is being considered as a differential diagnosis for insulin resistance in a diabetic, an investigation for acromegally (by IGF-1) may prove more fruitful.

Feline ACTH Stimulation Test

As with any test for hyperadrenocorticism (see canine HAC) this test may be affected by a variety of chronic illnesses causing a significant false positive rate.

- 1. Collect basal blood sample.
- 2. Inject 0.125 mg* or 5ug/kg of synthetic ACTH (Synacthen/Tetracosactide) i/v.
- 3. Collect two further blood samples at one and three hours later.
- 4 Label sample tubes clearly and request cortisol.
 - * 0.25mg may be used in cats over 5kg

Interpretation

Normal cats will increase up to about 400 nmol/L.

Feline Combined High-Dose Dexamethasone/ACTH Test

Because more than one diagnostic endpoint is assessed, this test may be more helpful than an ACTH stimulation test alone.

- 1. Collect basal blood sample.
- 2. Inject 0.1mg/kg Dexamethasone i/v.
- 3. Collect a second blood sample at 2 hours.
- 4. Immediately inject 0.125 mg of synthetic ACTH (Synacthen/Tetracosactide) i/v.
- 5. Collect a third blood sample at 3 hours (1 hour after ACTH).
- 6. Label sample tubes clearly and request cortisol.

Interpretation

Normal cats show at least 50% suppression to a value <40nmol/L after dexamethasone and normal cortisol response after ACTH stimulation (up to 400 nmol/L). Cats with HAC show little suppression after dexamethasone and an exaggerated response after ACTH.

ACTH

Endogenous ACTH is useful in differentiating PDH and ADH in cats but it should not be used as a diagnostic test in the first instance. Special sample handling procedures are required (see canine section, page 23).

Cortisol/Creatinine Ratio (CCR)

This is a very sensitive test to exclude HAC but must not be used to diagnose HAC, as it is very poorly specific because non-adrenal illness will commonly cause a positive result.

A morning urine sample is collected and this reflects the cortisol release over several hours. Samples taken from a litter tray in the home environment will have a lower false positive rate than those obtained from animals in the clinic. Non-absorptive material such as gravel, glass beads etc should be used as litter.

Normal cats have a CCR less than 10 x 10⁻⁶.

Cortisol/Creatinine Ratio (CCR) with Dexamethasone Suppression Test

When Cushing's disease is strongly suspected, this is a test which is especially useful in cats, as they usually get very stressed when visiting the vet. The owner does all the sampling at home and samples are sent to the lab together.

The procedure is the same as for the dog CCR/Dex. Test (page 24), except that the litter in the litter tray is replaced with washed and dried aquarium gravel and the urine samples collected with a plastic pipette. Wash the tray and gravel thoroughly between samples.

Guinea Pig Hyperadrenocorticism (HAC - Cushing's Disease)

Hyperadrenocorticism has been reported in Guinea Pigs. Affected animals display non-pruritic hair loss, thin skin, polydipsia, weight loss and decreased activity. The investigation of suspected cases is now much easier because of the availability of salivary cortisol measurement at NationWide Specialist Laboratories. Saliva samples can be collected on several large plain cotton buds before and after intramuscular injection of ACTH. The saliva cortisol approach eliminates the difficulties of multiple blood sampling in this small species and avoids unnecessary stress for the animal. See below for suggestions on collecting the saliva samples.

** Please note this protocol replaces the procedure in the Specialist Laboratory Services Information Manual, 19th Edition, March 2017, page 31.

Guinea Pig ACTH Stimulation Test

- 1. Collect basal saliva sample.
- Inject 5-10μg of synthetic ACTH IM
 (For a 1kg guinea pig: 0.04 mL Synacthen (250ug/ml) or 0.1 mL Tetracosactide (100ug/ml))
- 3. Collect 2 further saliva samples one and two hours later.
- 4. Place the saliva swabs in separate tubes and label clearly.
- 5. Request saliva cortisol analysis

Interpretation

Normal guinea pigs will increase up to about 725** nmol/L.

** Provisional range taken from Vet Record, June 23, 2007 p878-880

Notes on taking saliva samples.

Remove all food or eating material at least 2 hours before starting the test and taking the samples.

Use cotton wool buds to take the samples.

Use a carrot or similar tempting treat to generate some saliva production but do not give the treat until the last saliva sample is taken.

Ensure the cotton wool buds are really "wet" with saliva and take 2-3 if possible for each sample point. Place buds in a plain tube and label.

When the final sample is taken the Guinea Pig can have the tasty treat!

Equine Pituitary Pars Intermedia Dysfunction (PPID)

Equine Cushing's Disease (ECD) or Equine Cushing's Syndrome (ECS)

PPID is commonly known as equine cushing's disease/syndrome (ECD/ECS) and is caused by hyperplasia of the pituitary pars intermedia. This produces excess secretions of a variety of closely related peptides derived from pro-opiomelanocortin including alpha-MSH, B-endorphin and CLIP, and smaller amounts of ACTH, which results in increased adrenocortical production of cortisol. It is a disease of the older horse, usually >15 years. The most obvious clinical sign is the characteristic long curly coat (hirsutism) and abnormal shedding patterns. Other common clinical signs include weight loss, weight distribution (pot belly, loss of epaxial musculature, increased supraorbital fat) lethargy, laminitis, polyuria, polydipsia and hyperhidrosis. Single basal cortisol levels are often within or below the normal range unless the animal is acutely stressed due to pain possibly caused by an attack of laminitis. Insulin levels are frequently raised but insulin alone should not be used to diagnose ECS. Endogenous ACTH is considered to be a very sensitive test and can be done alone or in conjunction with either the regular overnight Dexamethasone test or the combined TRH/ overnight Dexamethasone test.

Overnight Dexamethasone Suppression Test

- 1. Take a basal blood sample at about 17.00 hours (5 pm).
- 2. Inject 0.04 mg/kg Dexamethasone i/m.
- 3. Take a further blood sample 20 hours later.
- 4. Label samples clearly and request cortisol.

Interpretation

Normal horses should suppress cortisol levels to below 30 nmol/L and ideally <20 nmol/L. Horses with high baseline values (>150 nmol/L) usually have less suppression and values below 40nmol/L would be considered normal.

Important Note: Normal horses may give false positive results if the test is done in the autumn months (August, September, October).

<u>Insulin</u>

The analysis of Insulin on a basal serum sample is very useful as this detects the presence of peripheral insulin resistance. Insulin levels are often very high (>250 μ U/mL) in ECS. Insulin levels in equine metabolic syndrome (EMS) are raised but there can be considerable variation over a 24-hour period.

Normal horses and ponies have insulin levels less than 65 µU/mL.

ACTH

In horses and ponies a single sample analysed for ACTH has been shown to be a very sensitive (100%) test for equine ECS. Normal ACTH levels vary according to the time of year and are significantly increased in the autumn (August, September, October) so this must be taken into account when interpreting results for the diagnosis of ECS. However very high levels (>300pg/mL) are still highly suggestive of ECS and a normal level in the autumn (<100 pg/mL) is more likely to be truly negative. Horses with EMS may have elevated ACTH due to the stress of their laminitic condition but levels will still not be as high as expected with ECS.

Cortisol/Creatinine Ratio (CCR)

This can be useful as a screening test to exclude ECS but must not be used to diagnose ECS as it is not very specific and other factors such as exercise may give a positive result. A morning urine sample is collected for analysis.

Normal horses have a CCR less than 20 x 10⁻⁶.

Equine Metabolic Syndrome (EMS)

Equine Metabolic Syndrome (EMS)

EMS was often previously referred to as Atypical Equine Cushing's disease or Peripheral Cushing's disease. EMS is a syndrome of obesity, insulin resistance (IR) and laminitis. Obesity is not a significant factor on its own and appears to be strongly linked to genetic predisposition especially in some of the susceptible British breeds of ponies. Obesity includes abnormal weight distribution including "cresty" neck or bulging supraorbital fat pads. Insulin resistance is the crucial factor in EMS and results in significant hyperinsulinaemia. Laminitis is the most common disease expression associated with EMS and is often seasonal, occurring during the summer months when animals have access to lush summer pastures. The spectrum of EMS varies from overly obese horses that become insulin resistant during the summer months resulting in seasonal hyperinsulinaemia and increased risk of summer pasture associated laminitis to the other end of the spectrum where insulin resistance and hyperinsulinaemia are present all year round and laminitis can occur at any time. Hirsutism or abnormal coat shedding are not usually found in EMS and this helps in the differentiation with Equine Cushing's Syndrome (ECS). Horses are also usually younger than those with ECS. Basal cortisol levels may be significantly increased above the reference range. Endogenous ACTH is expected to be within the normal range or borderline high. Animals test negative on an overnight Dexamethasone test, so this is a useful test to differentiate EMS and ECS except in the autumn when false positives may occur due to seasonal variation.

Diagnosis of EMS is obtained by confirming insulin resistance and showing hyperinsulinaemia in a basal serum Insulin sample ideally taken in the morning before feeding (fasted sample). Insulin levels are usually raised on basal samples, but if there is doubt, a dynamic test of insulin resistance can be used, e.g. the combined glucose-insulin test. Ensure horse is fasted or not fed any carbohydrate prior to collection of samples and the horse has not exercised that day. Elevated insulin values without having fed a carbohydrate meal is likely to be caused by insulin resistance.

In-feed Oral Glucose Challenge Test:

This test can be performed at the stable premises and is useful in suspected IR cases with normal fasted insulin concentrations:

- 1. Fast the patient overnight (12 hours).
- 2. Request owner to give a nonglycemic feed (e.g. chaff) containing 1g/kg glucose or dextrose powder (wet the feed to facilitata mixing and ingestion).
- 3. Blood sample for insulin 2 hours post consumption of the feed.

Interpretation

Serum insulin > 85 iu/l in the 2 hour sample is indicative of IR.

Adrenal Disease In Ferrets

Adrenal disease is commonly seen in ferrets over 3 years old and is caused by adrenal gland tumours often associated with the left adrenal gland (60% - 70%). Females appear to be more affected than males (70%). Clinical signs include significant hair loss, muscle wastage, aggression and vulval enlargement in females. Diagnosis can be confirmed by assaying the adrenal sex hormones, oestradiol (E2), 17-Hydroxy-Progesterone (OHP) (Ferret Adrenal Profile).

Additional diagnostic information may also be obtained by adding Androstenedione to the above hormones (Ferret Adrenal Profile PLUS). Cortisol alone is not a good diagnostic test for ferret adrenal neoplasia. Often the tumour has differentiated in such a way that cortisol is not its principal product.

Adrenal Sex Hormone Imbalance In Dogs (SHAP)

Adrenal sex hormone imbalance or congenital adrenal hyperplasia-like syndrome is a recently described disorder in dogs. The clinical signs are symmetrical alopecia and hyperpigmentation on the rump, perineum, caudal thighs, neck, tail and trunk, other areas being generally spared. Presentation is similar to the previously described growth hormone (GH) responsive alopecia and it is now thought that although decreased levels of GH may be a contributing factor, these are secondary to the sex hormone imbalance. The condition is thought to be due to a partial deficiency in glucocorticoid synthesis which results in the accumulation of the adrenal sex hormones that are precursors for these steroids. Unfortunately, such enzyme deficiencies have yet to be demonstrated in-vitro. Certain breeds are predisposed - Pomeranians and Chow Chows.

The pathophysiology of this presentation is somewhat controversial. Many of the previously described therapy responsive dermatopathies (GH responsive, castration response, Lysodren responsive, melatonin responsive etc) have very similar clinical presentations and perhaps their response to those therapies is non-specific. Realisation of these similarities has lead to the coining of the term "Alopecia X".

Adrenal sex hormones are also useful in assessing male dogs who are attractive to other dogs (male and female).

Recent studies have suggested that some dogs with hyperadrenocorticism (particularly those with adrenal neoplasia) do not test positive by conventional ACTH stim. test, but will show elevated levels of sex hormones, including OHP. It is suggested these dogs are also tested using the SHAP profile.

In some cases, SHAP testing may also be used to help identify stress responses which give false positive ACTH stimulation test results. The OHP levels may be less affected than the cortisol levels. Less than 10% of dogs with hyperadrenocorticism would be expected to have post ACTH OHP values <4.5 nmol/L (Chapman et al, Veterinary Record, December 2003).

SHAP Testing Protocol

Adrenal sex hormone status is investigated using a standard ACTH stimulation test (Page 22).

17-OH-Progesterone (OHP) has been identified as a sensitive marker and this along with Cortisol is measured before and after ACTH administration.

Interpretation

In normal dogs the basal OHP level is less than 1.0 nmol/L (usually undetectable (<0.1 nmol/L)) and after stimulation with ACTH normal dogs show an increase up to 2.0 nmol/L.

Dogs with a possible sex hormone imbalance often have a raised basal OHP level and show a significant increase after stimulation, greater than 3.5 nmol/L.

Dogs with HAC show a significant increase after stimulation. An exaggerated cortisol response could indicate classic hyperadrenocorticism.

Canine Full Adrenal Sex Hormone Profile

In addition to the 17-OH-Progesterone (OHP), we are now able to offer an expanded adrenal steroid hormone profile which includes Cortisol, OHP, Progesterone, E2, Testosterone and Androstenedione. This gives a complete picture of the adrenal status

A standard ACTH stimulation test (Page 22) should be performed. At least 2.0 mLs of serum or plasma is required to carry out all the analysis.

Aldosterone

Aldosterone is the most important mineralocorticoid and it is produced by the zona glomerulosa of the adrenal cortex. Release of aldosterone is mainly controlled by the reninangiotensin system (RAS) and by potassium levels in the blood.

Aldosterone levels may be useful in the differential diagnosis of primary and secondary Addison's and to assess the degree of destruction of the mineralocorticoid areas of the adrenals. An ACTH stimulation test should be performed and cortisol and aldosterone measured in both samples.

Hyperaldosteronism may be primary or secondary. Primary hyperaldosteronism (Conn's syndrome) is rare in the dog and cat and is usually caused by a small, solitary, aldosterone-producing adenoma of the adrenal cortex. Secondary hyperaldosteronism is more common and is caused by continued stimulation of the RAS.

Persistently low serum potassium levels <3.0 mmol/L (Hypokalaemia) may be due to hyperaldosteronism (aldosteronoma).

Because exogenous glucocorticoids will have little or no effect on the RAS, aldosterone production is not affected by them. By measuring aldosterone response to ACTH instead of (or in addition to) cortisol we can confirm primary hypoadrenocorticism in dogs that have already had several or more days of glucocorticoid therapy. Measuring only cortisol in the same circumstance risks a false positive diagnosis of hypoadrenocorticism because of the negative feedback effects of exogenous steroid therapy (including topical skin, eye and ear preparations) on endogenous glucocorticoid production.

Aldosterone levels can be measured in the same samples as cortisol and the measurement of aldosterone pre- and post- ACTH stimulation provides the most information on the mineralocorticoid producing abilities of the adrenals (see Hypoadrenocorticism: Diagnosis).

The principle uses of aldosterone analysis therefore are:

- Confirming primary hypoadrenocorticism in an animal which has already been receiving therapy prior to diagnosis
- Differentiating primary from secondary hypoadrenocorticism
- Investigating hypokalaemia (considering aldosteronoma)
- Assessing renin-angiotensin dependency in cardiac patients prior to ACE inhibition
- Assessing mineralocorticoid effect of anti-adrenal therapy (e.g. Lysodren, Vetoryl)

Example Aldosterone and Cortisol Interpretation:

	Cortisol (nmol/L)	Aldosterone (pmol/L)
Pre ACTH	<20	<20
Post ACTH	<20	<20

Consistent with primary hypoadrenocorticism (adrenocortical destruction)* – even if already received symptomatic glucocorticoid therapy

	Cortisol (nmol/L)	Aldosterone (pmol/L)
Pre ACTH	<20	Normal
Post ACTH	<20	Normal

Consistent with secondary (pituitary) hypoadrenocorticism OR

Exogenous glucocorticoids (including topical skin, eye and ear preparations)

^{*} In the case of trilostane (Vetoryl®) therapy, the ACTH response test results reflect the status of adrenal steroidogenesis only at the time of the test. Because of the temporary and reversible action of trilostane, the response may be greater later in the day, i.e., unlike in the case of mitotane, a flat aldosterone response does not necessarily mean that the trilostane patient is absolutely unable to make aldosterone.

Notes	®
	NationWide LABORATORIES
	01223 493400 www.thehormonelab.co.uk
a ⁱ	
Top -	199
KYM DYMAN	No. No.



Male Reproductive Function



Male Reproductive Function

Cryptorchid/Testicular Function

Testosterone

Testosterone may be analysed in a single sample if the presence of a whole testicle is suspected. Where testicular remnants or incomplete castration may have occurred it is better to do the hCG stimulation test (listed below).

hCG Stimulation Test

The hCG stimulation test may be used to identify cryptorchidism or testicular function in virtually any species. The following procedure may be used as a guideline; the product data sheet should be consulted for exact doses as required (*).

- 1. Take basal blood sample
- 2. Inject hCG i/v (Chorulon, Intervet) according to species:

Equine: 6000 IU

Canine: 200 - 500 IU*

Feline: 100 - 500 IU*

- 3. Take a second blood sample 30 minutes to 2 hours later.
- 4. Label both samples clearly and request testosterone.

Interpretation

Animals that have been castrated successfully have low basal levels of testosterone which show little or no increase in testosterone levels following the administration of hCG. Cryptorchid animals show increased basal levels of testosterone which usually increase following the administration of hCG. Animals with normal testicular function should have basal testosterone levels within the relevant normal range and stimulate adequately post injection of hCG.

Oestrone Sulphate (O-SO4) - Equine test for cryptorchidism (Rig Test).

Rig is a behavioural term used to describe a castrated animal showing unexpected male behaviour. A true rig is an animal whose unexpected male behaviour is confirmed to be due to testicular tissue. A false rig is an animal whose unexpected male behaviour is confirmed not due to testicular tissue. Oestrone Sulphate is a single sample test to distinguish between cryptorchids, true and false rigs. The O-SO4 is very accurate but any borderline results will need to be confirmed using the hCG stimulation test. O-SO4 must not be used in donkeys or animals under three years old.

Note: Both **O-SO₄ and the hCG stimulation test** should be used on any animals likely to be the subject of any dispute or court case.

AMH (Anti-Mullerian Hormone)

AMH may also be used as a cryptorchid test to detect testicular tissue and may also be used in young animals under 3 years old. AMH levels should be undetectable (AMH < 0.01ng/mL) in correctly castrated animals. Other tests (as above) should also be used on any animals likely to be the subject of any dispute or court case.

Notes	®
	NationWide LABORATORIES
	01223 493400 www.thehormonelab.co.uk
a ⁱ	
Top -	199
KYM DYMAN	No. No.



Female Reproductive Function



Female Reproductive Function

Ovulation Detection

Ovulation Detection in the Bitch

The use of quantitative blood progesterone estimations at specific times after the start of vulval bleeding can help to determine the best time for mating. This is very useful in bitches with reduced fertility or those who have "missed" on previous occasions due to unknown reasons.

Progesterone (quantitative)

A blood sample should be taken about day 7 post vulval bleeding. In most bitches the progesterone concentration will be <3 nmol/L. If the progesterone level is above 6 nmol/L, functional luteinisation has occurred and depending on the actual level we can advise when to sample again if required.

Interpretation

Progesterone <1.0 - 3.0 nmol/L Resample 3 - 4 days

Progesterone >6 nmol/L Indicates functional luteinisation, resample again if

required

Progesterone >20 nmol/L Ovulation has occurred, resample again if required

Progesterone 16 - 20 nmol/L Mate within 33 - 57 hours
Progesterone 20 - 38 nmol/L Mate within 9 - 33 hours
Progesterone >38 nmol/L Mate within 9 hours

Sample required – 0.5 mL serum (NO GEL TUBES) or heparin plasma.

Ovulation Detection in Other Species

Progesterone may be used to identify that ovulation has taken place in a wide variety of species by serial sampling to test progesterone concentration. Females in anoestrus will normally have undetectable progesterone levels. When ovulation occurs the progesterone rises and this rise can be detected with a sensitive quantitative progesterone.

Ovarian Tissue Detection

AMH (Anti-Mullerian Hormone)

AMH (Anti-Mullerian Hormone) is produced by ovarian granulosa cells in females and sertoli cells in males and so is a very useful diagnostic test to detect the presence of ovarian or testicular tissue.

Neutered status can be very difficult to determine, especially in animals with an unknown history such as strays or rescue animals. AMH is a very useful to determine spay status in both cats and dogs. Our studies have confirmed it is more reliable than performing the GnRH stimulation test where no signs of oestrus are present.

AMH levels remain constant throughout the cycle so females do not need to be showing signs of oestrus.

At NationWide Specialist Laboratories we have fully tested our "super-sensitive" AMH assay in conjunction with our local animal rescue centre, helping them to quickly assess the reproductive status of animals coming in to be re-homed

Our "super-sensitive" assay allows us to detect very low concentrations of AMH.

Animals that have been neutered will have undetectable levels of AMH (< 0.01 ng/mL).

Females who may have ovarian tissue present are likely to have detectable AMH levels greater than 0.01. However, negative results have sometimes been obtained from animals with a small remnant of ovarian tissue. The most common situation in which false negatives arise is during or shortly after (sometimes subclinical) oestrus. The additional measurement of progesterone (or oestradiol) will reduce the likelihood of missing remnant ovarian tissue when it is present.

A Progesterone assay should be also performed in an animal showing signs of oestrus where the AMH result was < 0.01 if there is still a strong suggestion of an ovarian remnant due to overt clinical signs.

0.5mL of serum is required.

Animals showing signs of oestrus at the time of sampling should have AMH and Oestradiol measured.

Animals who have been showing signs of oestrus 2-6 weeks previously at the time of sampling should have AMH and Progesterone measured.

Oestradiol (E₂)

Oestradiol may be used to detect the onset of oestrus in females or to detect ovarian tissue in spayed females. A single sample may only be required if the animal is showing signs of oestrus at the time of sampling. In spayed animals it may be necessary to do a GnRH or hCG stimulation test to stimulate any ovarian tissue to produce E₂. Males with suspected sertoli cell tumours should be screened with E₂.

Female Reproductive Function

GnRH Stimulation Test

- 1. Take a basal blood sample.
- Inject i/v 0.32 μg Buserelin (Receptal®, Intervet SP).
 (A very small volume of drug required (80μL) therefore it is advisable to dilute the drug.)
 - 1:10 with sterile water for injection and then inject 0.8mL).
- 3. Take a further blood sample 3 hours later.
- 4. Label tubes clearly and request Oestradiol.

Interpretation

Males and females in anoestrus should have undetectable (<10.0 pmol/L) levels of E2.

Ovarian activity or abnormal testicular activity will result in increased levels of E2.

Sertoli cell tumours in male dogs normally produce significant levels of oestrogens.

Each case must be carefully assessed and interpretation given according to clinical history.

Canine hCG Stimulation Test

- 1. Take a basal blood sample.
- 2. Inject i/v 200 500 IU hCG (Chorulon®; Intervet SP).
- 3. Take a further blood sample at 90 120 minutes
- 4. Label tubes clearly and request Oestradiol.

Interpretation

Basal oestradiol levels should be <10.0 pmol/L.

Levels greater than 10.0 pmol post hCG are suggestive of ovarian remnants.

Feline hCG Stimulation Test

This test is reported as being more effective at detecting ovarian remnants in queens. Testing should be done between one and three days after the onset of oestrous behaviour. The administration of hCG is to induce ovulation.

- 1. Take a basal blood sample (do not use serum gel tubes).
- 2. Inject i/m 500 IU hCG (Chorulon®; Intervet SP).
- 3. 7 days later take a further blood sample (do not use serum gel tubes).
- 4. Label tubes clearly and request Progesterone.

Interpretation

Basal progesterone levels should be <3.0 nmol/L.

After 7 days progesterone levels increase significantly above 3.0 nmol/L and often to >15.0 nmol/L in queens with ovarian remnants. Luteal tissue may be easier to see at exploratory surgery

Equine Inhibin

Inhibin is a polypeptide hormone secreted by ovarian granulosa cells in females and testicular sertoli cells in males. It selectively suppresses the secretion of pituitary FSH and also has local paracrine actions in the gonads.

The most common ovarian tumour in the mare is the granulosa cell tumour (GCT). Clinical diagnosis of GCTs is usually based on behavioural abnormalities such as prolonged anoestrus, stallion like behaviour, persistent oestrus and nymphomania and are often difficult to handle. Unilateral ovarian enlargement with atrophy of the other ovary is also a common finding.

Testosterone is elevated in approximately 54% of mares with GCTs and usually only those showing stallion-like behaviour. Mares showing anoestrous or persistent oestrous may have normal testosterone concentrations. Progesterone concentrations are almost invariably low (< 1.0 nmol/L) in mares with GCTs and as mares are usually acyclic, the detection of an elevated progesterone (>3.0 nmol/L) suggests that a GCT is unlikely to be present. Oestradiol levels are very variable and are not suitable as an aid to the detection of GCTs. Inhibin concentrations have been found to be elevated above normal values in approximately 87% of mares with GCTs and therefore appears to be a much more accurate indicator of the presence of a GCT than testosterone alone. The measurement of progesterone, testosterone and inhibin is suggested and we offer this as a GCT profile.

Equine AMH

Anti-Mullerian Hormone (AMH) is produced by granulosa cells in the equine ovary and in normal mares this level is very low (usually less than 1.0 ng/mL). AMH levels in normal mares do not significantly change during oestrus or pregnancy, unlike Inhibin, and this makes it clinically more useful in these cases.

Mares with GCTs have greatly increased levels of AMH (usually greater than 25ng/mL) and this makes it a very useful diagnostic tool in addition to Inhibin to detect the presence of GCTs. AMH analysis is especially important in pregnant or cycling mares where Inhibin and Testosterone levels are often raised as part of the normal cycle and in pregnancy.

We offer AMH analysis on its own, or as part of our GCT Profile which includes AMH and Progesterone.

Notes	®
	NationWide LABORATORIES
	01223 493400 www.thehormonelab.co.uk
a ⁱ	
Top -	199
KYM DYMAN	No. No.



Pregnancy Tests





Pregnancy Tests

Equine Progesterone(Early Indicative Pregnancy Test)

Quantitative progesterone is measured in a serum sample taken 18 - 25 days post service. Progesterone levels >13.0 nmol/L are highly suggestive of pregnancy at this time, levels <3.0 nmnol/L are likely to be non-pregnant. This is very useful for thoroughbred mares, where time is critical and who need to be returned to the stallion for further mating.

Progesterone can also be used in other ruminants that cycle every 21 days e.g. Llamas and Alpacas.

Note: Bovine samples for progesterone must be taken as clotted blood and separated as soon as possible. Heparinised bovine blood samples will have falsely low progesterone levels due to the red blood cells metabolising the progesterone.

Equine Pregnant Mare Serum Gonadotrophin (PMSG)

Note: SERUM sample MUST be used for PMSG

PMSG is produced from the endometrial cups from about day 40 - 120 post covering and is the next available pregnancy test. PMSG is prone to false positives as, if the foal dies, the endometrial cups continue to produce PMSG for some time. O-SO₄ should be used to confirm pregnancy and the presence of a live foal after 100 days.

Equine Oestrone Sulphate Pregnancy Test

Oestrone Sulphate can be used as a very sensitive pregnancy test in horses and demonstrates the presence of a live foal, if the foal dies the O-SO4 decreases very rapidly. $O-SO_4$ starts to rise from about 80 days post covering, increases steadily to above 25 ng/mL and remains detectable until very close to term.

 $O-SO_4$ may also be used in cows and donkeys as a late pregnancy test and in goats to distinguish between true pregnancy and "cloudburst". Pigs have an early $O-SO_4$ peak about 25-35 days and again at 80-110 days.

O-SO4 may also be used in alpacas who have a first peak at 21 – 28 days and also in the last trimester. Alpaca gestation is 343 days.

Total Urinary Oestrogens (TUE)

TUE may be used as a late pregnancy test (>100 days gestation) as there are very large concentrations of a variety of oestrogens excreted in the urine in late pregnancy. This is also applicable to a wide range of species.

Summary Table Of Equine Pregnancy Tests

Days Gest	Test	Use
18 - 25	Progesterone (quantitative)	Early - first indication
40 - 120	PMSG (quantitative)	Next - good indication
>100	Oestrone Sulphate	Middle/Late - very accurate
>100	Total Urinary Oestrogens	Middle/Late - urine test

Canine Pregnancy Test - Relaxin

Non-pregnant bitches have significant levels of progesterone as part of their normal cycle (metoestrus), therefore **progesterone alone cannot be used as a canine pregnancy test** but may be used to rule out pregnancy where very low levels are detected

Relaxin is very specific to pregnancy and so can be used as an accurate pregnancy test. Following conception relaxin rises to detectable levels from about 21 - 28 days. It is advisable to take blood samples at around 28 days or more post mating as this will allow for differences in conception dates.

Sample required is 0.5mL SERUM or HEPARIN plasma.

Sample must be collected and separated immediately and sent to lab on icepacks within 24 hours.

Feline Pregnancy Test - Relaxin

Progesterone may be analysed but this will only indicate if the cat is non-pregnant (Progesterone <3.0 nmol/L). Recent work has shown that feline relaxin may be used as a pregnancy test in the cat.

Faecal Pregnancy Detection

We have developed methods to determine pregnancy in wild cats such as lions, leopards and tigers using the measurement of progesterone in faecal samples. Please ask for details.

Notes	®
	NationWide LABORATORIES
	01223 493400 www.thehormonelab.co.uk
a ⁱ	
Top -	199
KYM DYMAN	No. No.



Gastrointestinal And Pancreatic



Gastrointestinal And Pancreatic

Trypsin-Like Immunoreactivity (TLI)

Serum TLI will detect exocrine pancreatic insufficiency (EPI) in dogs and cats in a single blood sample and this remains the test of choice for EPI in both species. Samples must be serum and the animals need to be fasted for at least 6 hours and preferably overnight. Acute pancreatitis gives rise to very high levels due to increased amounts of trypsinogen released into the general circulation. It is very important to sample suspected cases of pancreatitis as early as possible to obtain an accurate diagnosis, as the levels of TLI decrease as the condition progresses.

Interpretation

Animals with EPI have very low levels of TLI. Animals with pancreatitis may have very high levels of TLI.

Lipase (DGGR)

NationWide Specialist Laboratory are using a new improved method for pancreatic lipase activity which strongly correlates with cPLI in dogs and fPLI in cats.

In a study of 64 canine samples we derived cut-off values for the detection of pancreatic inflammation by comparison to Pancreatic Lipase Immunoreactivity (cPLI). This study confirmed excellent diagnostic agreement between "Lipase (DGGR)" and cPLI

- Overall agreement between "Lipase (DGGR)" and cPLI™ of 92%
- All positive cPLI cases (>400 ug/L; n=18) had "Lipase (DGGR)" > 95 IU/L
- All negative "Lipase (DGGR)" (<45 IU/L; n=30) cases had negative cPLI (<200 ug/L)
- 11 of 12 equivocal cPLI results (200 400 ug/L) were equivocal by "Lipase (DGGR)" (45 95 IU/L)
- Based on this study, cPLI appears to offer little diagnostic advantage over "Lipase (DGGR)"
- A similar smaller comparison study against fPLI yielded a cut-off of 25 IU/L

	Lipase (DGGR) (IU/L)	cPLI/fPLI (ug/L)
Canine normal	<45	<200
Canine probable pancreatic inflammation	>95	>400
Feline normal	<25	<3.6

 $cPLI = IDEXX Spec cPL^{TM}$ $fPLI = IDEXX Spec cPL^{TM}$

Vitamin B12/Folate

Serum vitamin B12 and folate concentrations are useful in diagnosing small intestinal malabsorption or "bacterial overgrowth" (now more commonly referred to as "antibiotic responsive diarrhoea" (ARD)) in dogs and cats. A single serum sample is required.

Interpretation

Low serum folate suggests malabsorption in the proximal small intestine. Low vitamin B12 suggests malabsorption in the distal small intestine. Low or normal vitamin B12 together with high folate suggests "bacterial overgrowth".

	Serum Folate (ng/mL)	Serum Vitamin B12 (pg/mL)
Decreased Level	Malabsorption/damage in the proximal small intestine.	Malabsorption/damage in the distal small intestine.
	Low dietary intake	Bacterial overgrowth (ARD)
	Drug interactions	Pancreatic insufficiency (EPI)
		Low dietary intake
		Hereditary selective malabsorption
Increased Level	Small intestine bacterial overgrowth (ARD) Pancreatic insufficiency (EPI)	High dietary intake
		Vitamin over supplementation
		Immunoproliferative disease
	High dietary intake	
	Vitamin over supplementation	

Pancreatic Lipase Immunoreactivity (PLI)

Canine and feline serum PLI (cPLI & fPLI) has the advantage of remaining raised long after the levels of TLI have decreased back into the reference range in pancreatitis cases.

Chronic small intestine disease which is complicated by chronic pancreatitis is more common in cats than primary chronic pancreatitis and so in these cases Vitamin B12 and Folate concentrations should also be measured.

Interpretation

Dogs: Normal: <200ug/L. Levels of cPLI >400ug/L are highly suggestive of pancreatitis. Cats: Normal: <3.5 ug/L. Levels of fPLI >5.4 ug/L are highly suggestive of pancreatitis

 $cPLI = IDEXX \ Spec \ cPL^{TM}$ $fPLI = IDEXX \ Spec \ fPL^{TM}$

Insulin

Insulin levels in normal fasting dogs are less than 29 μ U/mL. To investigate hypoglycaemia, samples for insulin assay should be taken during a period when the blood sugar is known to be low. If the insulin levels are greater than 29 μ U/mL this is considered to be excessive secretion. Guidelines for the interpretation of fasting insulin are given below.

Insulin greater than 29 µU/mL	High insulin, insulinoma very likely
Insulin 16 – 29 μU/mL	High/normal insulin, insulinoma likely
Insulin 8.3 - 16 0μU/mL	Low/normal insulin, insulinoma possible if glucose very low
Insulin less than 8.3 μU/mL	Undetectable/low insulin - insulinoma very unlikely

Repeat estimation of borderline results may be required to confirm the diagnosis of insulinoma. Elevated insulin in combination with normal or elevated glucose is consistent with insulin resistance rather than insulinoma.

Insulin/Glucose Ratio (IGR)

The IGR is useful when absolute hyperinsulinaemia (insulin >29 μ U/mL) is not present and the insulin levels alone are not diagnostic. Glucose and insulin are measured on the same sample.

$$IGR (Ulmol) = \frac{Serum Insulin (\mu U/mL)}{Serum Glucose (mmol/L)}$$

Interpretation

IGR of >7.0 U/mol is consistent with insulinoma when glucose is low or low-normal. An elevated ratio in an animal with high normal or high glucose indicates insulin insensitivity.

Amended Insulin/Glucose Ratio (AIGR)

The AIGR is reported to be more helpful in some cases and is calculated on the assumption that insulin levels should be undetectable if the blood glucose falls below 30 mg/dL (1.7 mmol/L).

AIGR (Ulmol) =
$$\frac{\text{Serum Insulin (}\mu\text{U/mL) }x100}{(\text{Serum Glucose (}mmol/L) }x18.02) - 30}$$

Note: If glucose is <30mg/dL then the divisor becomes 1.0.

Interpretation

AIGR of >50 is consistent with insulinoma when blood glucose is subnormal.

Equine Insulin

Horses with PIA and "peripheral hyperadrenocorticism" often have very high serum insulin levels due to peripheral insulin resistance.

Fructosamine

Fructosamine is useful to assess long-term glycaemic control in dogs and cats. It is a measure of plasma proteins which have undergone non-enzymatic glycation and is therefore related to the mean blood glucose concentration over the previous 1 - 3 weeks.

It is important that the fructosamine normal reference range is not used when monitoring diabetic patients as this is much lower. Results within the reference range usually indicate that the diabetic patient is likely to be experiencing significant periods of hypoglycaemia.

Haemolysed samples must not be submitted for the assay of fructosamine as this gives falsely elevated results.

Feline fructosamine levels greater than 550 µmol/L are consistent with poor diabetic control. Canine fructosamine levels greater than 450 µmol/L are consistent with poor diabetic control.

Notes	®
	NationWide LABORATORIES
	01223 493400 www.thehormonelab.co.uk
a ⁱ	
Top -	199
KYM DYMAN	No. No.



Calcium Disorders



Calcium Disorders


Parathyroid Hormone (PTH) And Ionised Calcium (iCA)

Parathyroid hormone (PTH) and ionised calcium (iCa) measurements are the mainstay of calcium investigations and their combined measurement should be considered in both hypoand hypercalcaemic situations where the cause is not immediately obvious. Neither total nor albumin corrected calcium measurements give a true reflection of calcium status in all cases. This is particularly so in animals with compromised renal function.

The combination of iCa and PTH measurement allow us to distinguish between four main categories of calcium regulation disorders:

- 1. Primary hyperparathyroidism (parathyroid dependent hypercalcaemia) functional parathyroid neoplasia
- 2. Parathyroid suppression (parathyroid independent hypercalcaemia) hypercalcaemia of malignancy, vitamin D toxicosis (cholecalciferol rodenticide, calcitriol, Dovonex), feline idiopathic hypercalcaemia* and granulomatous disease processes.
- 3. Primary hypoparathyroidism (parathyroid dependent hypocalcaemia) parathyroid gland destruction by inflammation or surgery
- 4. Secondary hyperparathyroidism (parathyroid independent hypocalcaemia) renal failure, calcium losses (e.g. pancreatitis, equine diarrhoea/colic), dietary deficiency, rickets, hyperadrenocorticism, hyperthyroidism.

*As it's name suggests, there is currently limited understanding of feline idiopathic hypercalcaemia. In some studies there is a suggested association with urinary acidifying diets. In others, calcium oxalate crystalluria has been reported and some amelioration achieved by a change to high fibre diets. Some authors believe the condition to be relatively benign, while others believe that the condition promotes the onset of renal dysfunction. Low-dose-glucocorticoids have been suggested as a therapy if there is concern about progressive disease.

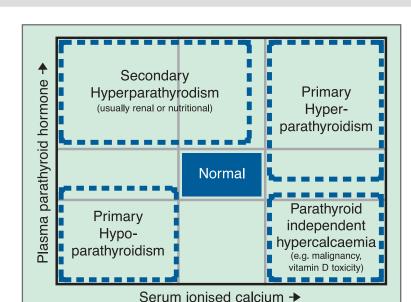

PTH/Ionised Calcium Interpretation Chart

The chart illustrates the diagnostic categorisation of calcium disorders using PTH and iCa alone. The chart as depicted applies well to dogs but the lower and smaller feline PTH reference range makes interpretation sometimes a little less clear-cut in cats.

Sample preparation is critical for the accurate analysis of PTH, see procedure below.

An ionised calcium sample should also be taken at the same time as the PTH sample as this gives an accurate indication of the calcium status at that time.

PTH Sample Preparation

Submitted sample must be **EDTA plasma or EDTA plasma with Aprotinin**.

- 1. Request a transport pack for delivery to the laboratory.
- Take the blood sample into a fridge cooled EDTA or Aprotinin EDTA tube (supplied in pack). Note: To avoid agglutination of samples in the assay, if possible, half fill several EDTA tubes with blood to ensure a high concentration of EDTA, especially if very high levels of calcium have been observed.
- 3. **Mix very well but gently at least 10 times** and centrifuge as quickly as possible. **Note:** If the Aprotinin EDTA tube will not fit into your centrifuge, then carefully transfer the well mixed blood into cooled tubes that do fit into the centrifuge.
- 4. Transfer the plasma into a cooled plastic (not glass!), EDTA tube which has been cooled in the fridge.
- 5. Immediately freeze (<-10°C) the plasma sample and keep frozen until dispatch in the transport pack. Freeze the plasma sample first, **separately**, before placing in the frozen pack for transport. Samples are stable for several weeks frozen.
- 6. Mark the tube with "EDTA plasma" and animal/owner name.

Ionised Calcium (iCa)

Submitted sample must be separated serum. EDTA plasma is not acceptable, the EDTA would chelate all the available calcium making it unavailable for analysis.

If possible take the sample at the same time as the PTH sample.

Collection technique

- 1. Obtain a blood sample in a plain **non-gel** tube and fill to the brim with blood.
- 2. Screw the lid down so no air enters the tube and allow to clot at room temperature for at least 30 minutes.
- 3. Spin the sample and separate the serum into another plain tube, take care to leave as small an air gap as possible.
- 4. Mark the tube with "Serum" and animal/owner name
- 5. The iCa separated serum sample may be frozen or stored in the fridge and send with the PTH sample outside the pack.
- 6. The laboratory will apply a correction formula to the result to take account of the changing pH that will have occurred due to the exposure to air. This will result in an estimated iCa at a standardised pH.

Parathyroid Hormone Related Peptide (PTHRP)

PTHrP may be used for the differential diagnosis of hypercalcaemia of unknown origin where other diagnostic tests have not identified the aetiology of the hypercalcaemia but where PTH and iCa suggest parathyroid independent hypercalcaemia. PTHrP is a hormone that can be produced by dogs, cats and horses with several different types of tumours and is considered to be the underlying cause of hypercalcaemia of malignancy in many but not all cases. Almost all anal gland apocrine neoplasia, the majority of hypercalcaemic lymphomas and smaller percentages of myeloma and carcinomas are PTHrP positive. Circulating levels of PTHrP in normal dogs are almost undetectable (< 0.1 pmol/L). Levels greater than 0.1 pmol/L are considered potentially significant in dogs. Cats appear to show similar values.

Please note: Sample preparation is <u>critical</u> for the accurate analysis of PTHrP (see procedure for PTH, page 61). The Aprotinin EDTA tube supplied with the transport pack <u>should</u> be used for the collection of PTHrP samples.

Vitamin D Metabolites

Two forms of vitamin D are measurable in veterinary serum samples: 25-hydroxyvitamin D (25OHD) and 1,25-dihydroxyvitamin D (calcitriol).

25OHD is produced by the liver and its concentration parallels that of available vitamin D that may be dietary in origin (cholecalciferol (D3), ergocalciferol (D2)) or produced in the skin of certain animals under the influence of UV light. Because the enzymatic hydroxylation of vitamin D in the liver depends almost entirely on the availability of the substrate, this test is an excellent marker of overall vitamin D status. It can be used to diagnose both conditions of vitamin D deficiency and excess and therefore is valuable in the investigation of both hyper- and hypocalcaemic disorders. 25OHD is also the analyte of choice in verifying adequacy of UV light exposure in reptilian species. Other situations in which we can expect abnormal 25OHD results would include dietary deficiencies, malabsorption syndromes and calciferol-rodenticide toxicities.

Calcitriol is produced by renal tubular cells as a result of enzymatic action (1α-hydroxylase) on 25OHD substrate. The rate of this process is controlled by PTH concentrations (increasing PTH causes increased 1α-hydroxylase activity). Calcitriol is the most biologically potent form of vitamin D and its principle activities are directed at increasing serum calcium concentrations, including increased intestinal uptake of calcium. The failure of renal tubular cells to generate calcitriol in renal disease is one of the contributing mechanisms to renal secondary hyperparathyroidism (and "rubber jaw"). The measurement of calcitriol may be of some value in understanding, renal tubular function and the effects of PTH on the vitamin D system.

In dogs, cats and most veterinary species, very little passive intestinal absorption of calcium occurs. In these species, intestinal calcium absorption is almost exclusively mediated by vitamin D. The situation is different in horses, rabbits (and the hippopotamus!), in which there is significantly more passive absorption of calcium and calcium status is controlled more by renal losses. In these species, the effects of renal disease on calcium status will be different than in most common mammalian species.

Recent research has suggested that low concentrations of 25OHD may be a risk factor for congestive heart failure (CHF) in dogs. There is extensive evidence that Vitamin D has cardioprotective actions by suppressing the renin-angiotensin-aldosterone system, decreasing myocardial hypertrophy, inhibiting pro-imflammatory cytokines and improving endothelial dysfunction and atherosclerosis. (Kraus, Rassnick, Wakshlag et al, J Vet Intern Med, 2014, 28, 109 – 115).

Animals should be fasted overnight.

25-Hydroxyvitamin D sample is min 0.5mL separated serum or plasma. Shipping on ice is preferred and avoid excessive exposure to sunlight.

Calcitriol sample is min 0.75 mL separated serum. Shipping on ice is preferred and avoid excessive exposure to sunlight.

Notes	®
	NationWide LABORATORIES
	01223 493400 www.thehormonelab.co.uk
a ⁱ	
Top -	and the same of th
KYM DYMAN	No. No.



Other Endocrine



Other Endocrine

Insulin-Like Growth Factor 1 (IGF-1)

IGF-1 is produced by the liver in response to circulating growth hormone concentrations. Growth hormone (GH) is released in a pulsatile manner, has a short half-life and requires species specific assays. IGF-1 on the other hand, is stable in plasma, has a long half-life and does not require species specific assays (just species validated assays). The properties of IGF-1 make it an excellent markers for overall (average) GH output as IGF-1 levels correlate very closely with GH secretion. This enables an estimation of the GH status to be made in dogs and cats. IGF-1 is very stable and no special sampling procedures are required.

IGF-1 may be used to diagnose conditions of growth hormone deficiency (pituitary dwarfism), growth hormone excess (pituitary or progestogen/metoestrus associated acromegaly), monitoring acromegaly therapy and in monitoring dietary sufficiency (objective documentation of malnourishment, neglect).

Other factors to consider when interpreting IGF-1 are that concentrations correlate positively with breed size in dogs and inversely with age. It can be helpful to have a sample from an age/breed matched control (e.g. littermate) particularly when investigating suspected pituitary dwarfism. This is especially important when sampling small or "dwarf" breeds of dogs and cats. The laboratory does have access to a limited list of IGF-1 results from normal animals of varying breeds and ages to aid interpretation on request.

The narrower range feline breed sizes makes this less of an interpretative problem for cats unless investigating "dwarf breeds". However, animals that are unwell may have depressed IGF-1 concentrations. This may be important when investigating suspected acromegaly in a poorly control diabetic cat. IGF-1 concentrations >1000, are consistent with acromegaly in cats. However, in a newly diagnosed cases, values >700 could be suspicious and warrant retesting after improved diabetic control, this is because of the suppressive effects of "illness" associated with recent diabetes diagnosis or poor control.

When documenting malnutrition, paired samples before and a few weeks after instituting an adequate diet are useful.

IGF-1 has no value in documenting "growth hormone responsive" dematopathies.

Erythropoietin (EPO)

Erythropoietin (EPO) plays a key role in the regulation of red cell mass and erythrocyte production. EPO may be used to aid the differential diagnosis of polycythemia and non-regenerative anaemia. Raised levels of EPO occur in secondary polycythemias and renal tumours. Low/normal levels may occur in renal failure and polycythemia vera.

Gastrin (Special Request Assay)

Gastrin measurement is useful in the investigation of chronic vomiting in dogs where gastrinoma is considered a differential diagnosis. High fasting serum concentrations or an exaggerated response to feeding can indicate Gastrinoma in the dog.

"False" elevation in Gastrin can occur with certain medications including proton-pump inhibitors (e.g. Omeprazole) and in chronic renal failure.

Sample Preparation for the Assay of Serum Gastrin

Fast the animals overnight or for 12 hours or more.

Take the blood into a plain tube with no additive and allow to clot and centrifuge and separate off the serum within 1 hour of collection.

The samples may be stored for up to 8 hours at 2-8°C or below -20°C for up to 2 weeks. Avoid freeze-thaw cycles and keep frozen until assayed.

Samples must be sent to the laboratory frozen and courier is essential.

Lipaemic samples may interfere with the assay.

At least 1.0 mL of Serum is required for the assay.

Renin (Special Request Assay)

Renin is released by the juxtaglomerular apparatus of the kidney in response to a decrease in arterial blood pressure (including decreased blood volume) or a reduced Na content of glomerular ultrafiltrate. It is the first part in the Renin-Angiotensin-Aldosterone-System

Renin activates Angiotensin by cleaving Angiotensinogen to release Angiotensin I which through the action of Angiotensin converting enzyme (ACE) produces Angiotensin II. In turn, Angiotensin II, stimulates the release of aldosterone from the adrenal cortex. The action of aldosterone should promote the resorption and retention of Na and restore blood volume.


The measurement of Renin can be particularly useful in the investigation of hypokalaemia and hypertension.

In conditions where there is an elevated aldosterone it is helpful to discover if that is "appropriate" and is the result of Renin activation or "inappropriate" and likely to be autonomous such as functional adrenal cortex neoplasia (aldosteronoma).

The measurement of Renin is therefore essential in the investigation of hypertension or

situations of elevated aldosterone.

We now offer a full Hypotension/Hypokalaemia profile comprising Plasma RENIN, Aldosterone and Na/K ratio including the calculation of Aldo/Renin ratio.

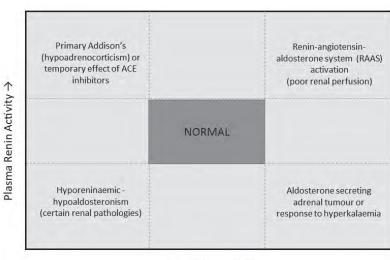

Aldosterone →

Sample preparation for the assay of plasma Renin


Blood must be collected into **pre-chilled** tubes containing EDTA as an anticoagulant. HEPARIN must not be used as this interferes with the assay.

The samples must be kept cold and ideally centrifuged in a refrigerated centrifuge to separate the plasma. The plasma samples must then be aliquoted and frozen at less than - 10°C until assayed.

Samples must be sent to the laboratory frozen and courier is essential.

Haemolysed or lipaemic samples must not be used in the assay so great care is needed when blood sampling to ensure samples do not become haemolysed.

At least 2.0 mL of EDTA plasma MUST be sent for the analysis of Renin otherwise it may not be possible to get a valid result.

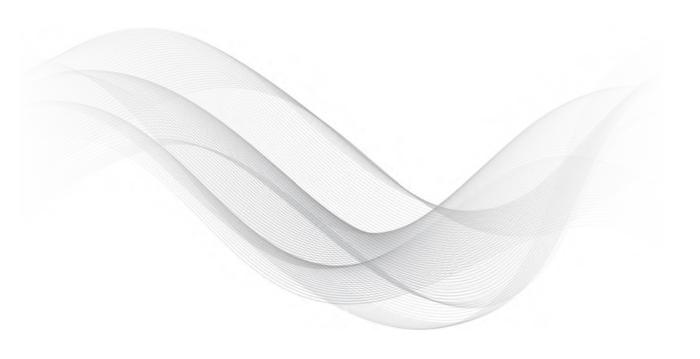

Glucagon (Special Request Assay)

The measurement of glucagon is helpful in the investigation of unusual cases of insulin resistant diabetes mellitus and in clinical research situations focused on factors affecting glucose metabolism.

Glucagonomas have been identified in dogs. A dermatological presentation has been associated with Glucagonoma in humans also known as necrolytic migratory erythema or hepato-cutaneous syndrome. In dogs, this skin condition usually has a hepatic origin.

Sample Preparation for the Assay of Plasma Glucagon

Blood must be collected into <u>iced</u> glass tubes containing EDTA as an anticoagulant. The tubes <u>must</u> be immersed in an ice-bath at 2 – 8°C throughout the entire collection and handling procedure.

Immediately following collection add 5000 Kallikrein Inactivating Units (KIU) per 10mL of whole blood in the form of Trayslol (500uL) or Aprotinin (1mg). Mix gently but thoroughly and keep chilled (on ice) until separation.

Centrifuge the samples in a refrigerated centrifuge at 2000g to separate the plasma. Separate the plasma using a glass pasteur pipette and transfer the plasma into glass tubes (aliquots of 0.6mL required) and freeze immediately. Store the plasma frozen at less than -10°C until assayed.

Samples must be transported to the laboratory frozen and remain frozen until assayed. Dry ice is recommended. Lipaemic and haemolysed samples may affect the assay so great care is needed when blood sampling to ensure samples do not become haemolysed. Courier transport is essential.

At least 1mL of EDTA plasma is required for the assay.

PLEASE NOTE AT THE CURRENT TIME WE ARE UNABLE TO PERFORM GLUCAGON ANALYSIS DUE TO THE WITHDRAWAL OF THE METHOD WE HAD VALIDATED.

WE WILL BE LOOKING AT OTHER METHODS FOR THE FUTURE.

NT-ProBNP - Canine/Feline Cardiac Marker

Heart disease in dogs and cats may be congenital (present at birth) or acquired. Acquired heart disease is more common and tends to occur as the animal gets older. Often the owner does not notice as severe clinical signs may not arise during gradual deterioration of cardiac function.

NT-pro-BNP will be a useful screen where subtle clinical signs are present, as part of a geriatric profile or a health check for breeds who have a history of cardiac problems.

Canine Interpretation

NT-pro-BNP levels less than 900 pmol/L are found in normal healthy dogs and it is unlikely that the clinical signs are related to cardiac failure.

NT-pro-BNP levels between 900 and 1800 pmol/L are in the suspect range and in the absence of a murmur or common clinical signs, heart disease is unlikely.

NT-pro-BNP levels greater than 1800 pmol/L Heart disease is likely. Further cardiac workup or specialist cardiac referral suggested.

Feline Interpretation

NT-pro-BNP levels less than 100 pmol/L are found in normal healthy cats and it is unlikely that the clinical signs are related to cardiac failure.

NT-pro-BNP levels between 100 and 270 pmol/L are in the suspect range and in the absence of a murmur or common clinical signs, heart disease is unlikely.

NT-pro-BNP levels greater than 270 pmol/L Heart disease is likely. Further cardiac workup or specialist cardiac referral suggested.

In cats and dogs NT-pro-BNP levels may decrease with cardiac medication. Arrhythmias and the presence of pulmonary hypertension may result in higher NT-pro-BNP levels. Azotemic animals may have increased NT-pro-BNP levels.

Please note very low values are suggestive of inappropriate sampling or sample degradation and a repeat sample should be submitted.

NT-pro-BNP Sample Preparation

Sample must be at least 0.5mL of separated EDTA PLASMA
SAMPLES MUST BE SEPARATED FROM THE RED BLOOD CELLS WITHIN 10
MINUTES.

Send to lab immediately (samples are stable for 24 hours at room temperature).

If there is going to be any delay then send to the lab on ice packs.

Do not send lipaemic or haemolysed samples.

Canine Benign Prostate Hyperplasia (BPH)

BPH appears to affect over 50% of entire male dogs over 5 years old. There is an increase in prostate cell size and cell numbers. Clinical signs suggestive of BPH are constipation, dysuria, haematuria, pain, lameness, behavioural problems, infertility and haematospermia. A simple blood test is now available and can provide a reliable alternative to rectal palpation. A single serum sample overcomes the many pitfalls of monitoring by rectal palpation such as the unpleasantness for dog, owner and vet, the awkwardness in large and miniature breeds, difficulty in distinguishing normal from abnormal and the fact that BPH must be severe to detect accurately by palpation.

Canine Prostate Specific Esterase (CPSE) is a major protein secreted by the prostate and represents more than 90% of the proteins found in seminal fluid. CPSE passes into the bloodstream following either an increase in its production or an alteration in the blood-glandular barrier caused by changes in glandular structure.

Interpretation

Normal: < 50 ng/mL. Strong suspicion of BPH: > 70 ng/mL

BPH Sample Preparation

Sample must be SERUM

Normetanephrine/Creatinine Ratio for Phaeochromocytoma

Diagnosis of phaeochromocytoma (PC) is based on a combination of clinical suspicion, finding an adrenal mass and measuring plasma or urine catecholamine metabolites such as normetanephrine. Confirmation is by histopathology. Affected dogs may present with weakness, polyuria/polydipsia, panting, restlessness, anxiety, weight loss or muscle wasting. Consequently, hyperadrenocorticism is a common differential diagnosis and finding an adrenal mass can cause confusion without further laboratory tests.

In canine patients, urine normetanephrine can be easily measured in an acidified urine sample.

Normetanephrine/Creatinine ratio (NCR) is the best parameter to differentiate dogs with PC from healthy dogs, dogs with hyperadrenocorticism or non-adrenal disease.

Dogs with PC have a very high NCR (> 400), compared to normal young dogs (< 100) or those with hyperadrenocorticism or non-adrenal disease.

Sample requirement: 10mL of acidified urine.



Therapeutic Drug Monitoring



Therapeutic Drug Monitoring

Phenobarbitone

Phenobarbitone should be given for about 2 - 3 weeks to allow for serum levels to stabilise before checking therapeutic concentrations.

Monitoring of therapeutic levels of phenobarbitone should be assessed at 8 - 12 hours post treatment. Peak levels occur about 5 hours post treatment.

Therapeutic levels for phenobarbitone are 15 - 40 µg/mL

SERUM GEL TUBES MUST NOT BE USED FOR PHENOBARBITONE ANALYSIS.

Potassium Bromide

Serum levels of potassium bromide may take at least 125 days to stabilise due to the long half life of the drug. Monitoring of therapeutic levels of potassium bromide can be assessed at any time post dosing.

Therapeutic levels for potassium bromide are 0.5 - 1.9 mg/mL. Young healthy dogs may tolerate levels up to 2.5 mg/mL and most dogs will tolerate levels up to 1.5 mg/mL.

SERUM GEL TUBES MUST NOT BE USED. EXTRA CARE IS NEEDED TO AVOID HAEMOLYSIS AS THIS SEVERLY ELEVATES THE RESULTS. SAMPLES MUST BE CENTRIFUGED AND SEPARATED BEFORE DESPATCH. WE SUGGEST THE USE OF LITHIUM HEPARIN SAMPLES AS THIS ALLOWS IMMEDIATE CENTRIFUGATION.

Potassium Bromide samples are reported with a haemolysis index where appropriate.

Levetiracetam (Keppra®)

Levetiracetam can be administered in addition to potassium bromide and phenobarbitone to further help in the reduction of seizure frequency and duration. The use of Levetiracetam is not associated with additional side effects and may help in cases of liver toxicity as the dose of KBr may be reduced as part of combination therapy.

A therapeutic range has not been established for the dog or cat but the human therapeutic range for Levetiracetam of 29 – 123 umol/L is considered appropriate. A steady state is not achieved with Levetiracetam, because accumulation does not occur. However the pharmacodynamic effect is believed to outlive the known half life. The risk of toxicity is believed to be low.

To establish half-life in an individual animal, sample at both peak (2-3 hrs post pill) and trough (immediately before next scheduled dose). Estimation of half-life is particularly relevant when being used along with phenobarbital because of its effect on levetiracetam metabolism.

Single sample monitoring should be trough (just before the next scheduled dose).

SERUM GEL TUBES MUST NOT BE USED FOR LEVETIRACETAM ANALYSIS.

Digoxin

Digoxin should be given for about 3 - 5 days to allow for serum levels to stabilise before checking therapeutic concentrations.

Monitoring of therapeutic levels of digoxin should be assessed at 8 hours post treatment.

Digoxin concentrations of 0.8 - 2.0 ng/mL are required post treatment.

SERUM GEL TUBES MUST NOT BE USED FOR DIGOXIN ANALYSIS.

Cyclosporine

Cyclosporine has been shown to be useful in the treatment of anal furunculosis (AF) and atopy in dogs.

Monitoring of cyclosporine levels during treatment may be useful. Approx 30 days should be allowed to reach stable trough levels using combination therapy. Approx 17 days should be allowed to reach stable trough levels using Cyclosporine alone.

Therapeutic levels of Cyclosporine about 12 hours post pill are expected to be between 300–400 ug/L for AF but lower for atopy.


Sample requirement is at least 2.0 mL of whole EDTA blood.

Zonisamide

Zonisamide has been used on its own and in combination therapy with Potassium Bromide. Routine monitoring of Zonisamide is reccommended for efficiacy and safety.

The human reference range of 10 - 40 ug/mL is commonly used in dogs.

SERUM GEL TUBES MUST NOT BE USED FOR ZONISAMIDE ANALYSIS.

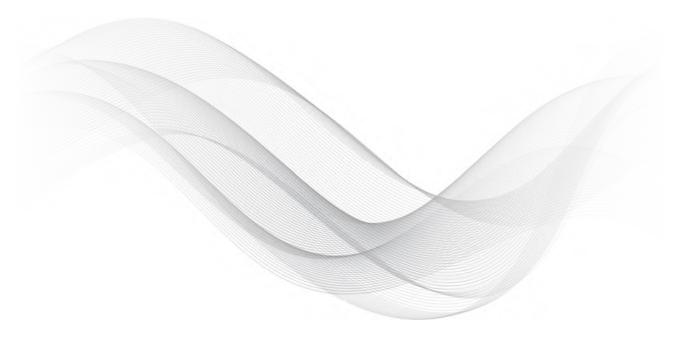

Notes	
	- NationWide
	LABORATORIES 01223 493400 www.thehormonelab.co.uk
	www.tileiloffiloflelab.co.uk
	N ROSE
1000 B x 17-18 9 0	
AND THE RESERVE OF THE PARTY OF	JUL MO - M



Endocrine Therapeutic Monitoring



Endocrine Therapeutic Monitoring

Basic Vetoryl® Monitor

ACTH Stimulation and Na:K ratio

In order to monitor Vetoryl® it is important to know what effect the medication is having on adrenal function. The test is performed 4-6 hours after the medication when Trilostane is at its greatest level of activity. By sampling at this time we know what the greatest level of adrenal suppression will be. It is likely that later in the day, suppression will be less. There are circumstances when investigating problem cases, that it may be appropriate to perform the ACTH stimulation test outside this period on the guidance of a specialist or Dechra Technical Support (01939 211200).

Vetoryl® will inhibit both glucocorticoid and mineralocorticoid production. Sodium and potassium are included because they will help identify whether there is too much suppression of the adrenal glands.

In interpreting the ACTH stimulation test 4-6 hours post Vetoryl® Capsule, we are looking for evidence of suppressed steroidogenesis, that is, post-ACTH cortisol concentrations <250nmol/L that also show an increase over baseline levels. The failure of ACTH to stimulate an increase in cortisol concentrations suggests that consideration be given to lowering the dose or withholding medication.

Basic 'Pre Pill' Vetoryl® Monitor

Single Cortisol and Na:K ratio

Vetoryl® Monitor Plus

ACTH Stimulation, Na:K ratio, Alkaline Phosphatase (AP), Alanine Aminotransferase (ALT), Urea, Creatinine, Phosphorus, Total Protein, Albumin, Globulin, Cholesterol and Glucose

The Vetoryl® datasheet recommends regular monitoring particularly during the early phase of therapy that includes ACTH stimulation and biochemistry including electrolytes.

The Vetoryl® Monitor Plus includes the ACTH stimulation and Sodium:Potassium ratio but adds a tailored selection of other biochemistry tests.

'Pre Pill' Vetoryl® Monitor Plus

Single Cortisol, Na:K ratio, Alkaline Phosphatase (AP), Alanine Aminotransferase (ALT), Urea, Creatinine, Phosphorus, Total Protein, Albumin, Globulin, Cholesterol and Glucose

	,
Alkaline phosphatase (AP) Alanine aminotransferase (ALT)	Both AP and ALT are usually elevated in untreated hyperadrenocorticism. AP is elevated in dogs because of induction of a steroid induced isoenzyme and ALT is elevated because of a "steroid hepatopathy". By measuring these enzymes, it is possible to assess improvements in liver pathology and circulating steroid concentrations. The Vetoryl® datasheet recommends that dogs should be monitored at regular intervals for primary hepatic disease.
Urea Creatinine Phosphorus	The measurement of urea, creatinine and phosphorus will assist in the identification of renal disease (and iatrogenic hypoadrenocorticism). Renal dysfunction is not uncommon in geriatric dogs. The Vetoryl® datasheet recommends that dogs should be monitored at regular intervals for renal disease and that the product should not be used in dogs with renal insufficiency.
Total protein Albumin Globulin (calculated)	The measurement of proteins aids in the assessment of hydration status in dogs that may be polydipsic. Albumin and globulin measurement can help identify late stage (protein-losing) glomerulopathies and inflammatory disorders that can affect dogs with hyperadrenocorticism.
Cholesterol	Cholesterol is commonly elevated in dogs with untreated hyperadrenocorticism. Measurement during therapeutic monitoring provides information on whether deranged lipid metabolism has been corrected.
Glucose	Diabetes mellitus is not uncommon in dogs with both untreated and treated hyperadrenocorticism. The Vetoryl datasheet recommends that dogs are monitored at regular intervals for diabetes mellitus.



Basic Thyforon®

Thyroxine (T4)

The measurement of thyroxine (T4) at a single point in time. It is important to know that a dog receiving Thyforon® therapy is absorbing the medication appropriately. Insufficient absorption or insufficient dose will lead to treatment failure and over-dose can lead to clinical signs of hyperthyroidism. When a T4 measurement is made, it is important to record the hours post-pill at which the sample was taken to be able to appropriately interpret the results.

The above monitoring may also be applied for other canine thyroid medication such as Leventa® and Thyroxanil®

Thyforon® Monitor Plus

Thyroxine (T4) and Thyrotropin (TSH)

The measurement of thyroxine (T4) and thyrotropin (TSH) at a single point in time. It is important to know that a dog receiving Thyforon® therapy is absorbing the medication appropriately. Insufficient absorption or insufficient dose will lead to treatment failure and over-dose can lead to clinical signs of hyperthyroidism. When a T4 measurement is made, it is important to record the hours post-pill at which the sample was taken to be able to appropriately interpret the results.

Thyrotropin (TSH) levels take several days to change in response to changes in thyroxine (T4) therapy, they therefore are able to reflect the adequacy of therapy in the days preceding the test. TSH could be described as the "fructosamine of hypothyroidism".

The combined measurement of T4 and TSH will help in the identification of compliance failure and reduce unnecessary dose adjustments that could arise from reliance on T4 measurement alone on an "unrepresentative" day.

The above monitoring may also be applied for other canine thyroid medication such as Leventa[®] and Thyroxanil[®].

Basic Felimazole®

Thyroxine (T4)

The measurement of single thyroxine (T4) concentration in hyperthyroid cats allows us to assess the adequacy of therapy and to make appropriate dose adjustments. As long as the cat is receiving medication on its regular schedule the interval post-pill at which the sample is taken is not important.

The above monitoring may also be applied for other feline hyperthyroid medication such as Vidalta[®] (MSD), Thyronorm[®] (Norbrook), Thiafeline[®] (Animalcare), Tapazole[®] (King Pharmaceuticals).


Felimazole® Monitor Plus

Thyroxine (T4), Alkaline Phosphatase (AP), Alanine aminotransferase (ALT), Urea, Creatinine and Phosphorus

Felimazole® monitor plus includes the measurement of thyroxine (T4) but also includes specific biochemistry tests.

Alkaline phosphatase (AP) Alanine aminotransferase (ALT)	Both AP and ALT are commonly elevated in untreated hyperthyroidism. By measuring these enzymes, it is possible to assess improvements due to therapy and to identify any underlying primary hepatic disease.
Urea Creatinine Phosphorus	The measurement of urea, creatinine and phosphorus will assist in the identification of renal disease. The "unmasking" of previously subclinical renal insufficiency is not uncommon following all forms of anti-thyroid therapy. By monitoring these parameters, it is possible to better navigate the balance between hyperthyroidism and renal insufficiency in case this situation should arise.

The above monitoring may also be applied for other feline hyperthyroid medication such as Vidalta® (MSD), Thyronorm® (Norbrook), Thiafeline® (Animalcare), Tapazole® (King Pharmaceuticals).

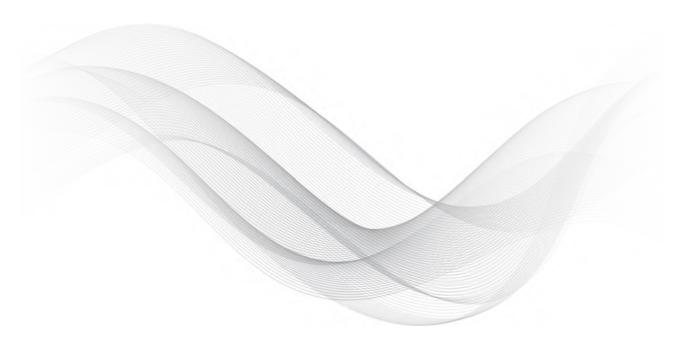





Serology/ Immunology



Serology/Immunology

Acetylcholine Receptor Antibodies (ACRAB)

Acquired canine myasthenia gravis (cMG) is an autoimmune disease of neuromuscular transmission caused by autoantibodies directed against the acetylcholine receptors at the neuromuscular junction. cMG can occur as two forms, focal and generalised. In focal MG dogs usually present with megaoesophagus as weakness affects the oesophageal, pharyngeal and facial muscles. In generalised cMG there is widespread skeletal muscle weakness which is made worse by exercise. ACRABs may only be used to diagnose the acquired form of the disease.

Positive serum titre for ACRABs is highly suggestive of acquired cMG.

A 1mL serum sample is required.

Type 2M Muscle Fibre Antibody (2MAB)

Disorders of the muscles of mastication occur relatively commonly in clinical practice and can be of myopathic or neuropathic origin. Masticatory muscles are predominantly composed of fibres designated as Type 2M. Masticatory muscle disorders (MMD) can be diagnosed by the assay of serum from affected animals to see if there are autoantibodies against Type 2M fibres.

False negative results may be obtained if immunosuppressive dosages of corticosteroids have been given for longer than 7 - 10 days in end-stage masticatory muscle myositis and in polymyositis. A muscle biopsy would be necessary for confirmation of the diagnosis in these cases.

A 1mL serum sample is required.

Sarcoptic Mange Antibody (SMA)

It is useful to be able to make a firm diagnosis for the presence of sarcoptic mange mites as this helps to choose the right treatment regime especially where there may be breed sensitivity to some drugs. Skin scrapes often fail to detect the presence of mites due to the difficulty in actually finding the mites.

A serum sample taken 4 - 8 weeks post infection will detect the presence of sarcoptic mange antibodies (IgG). A borderline result should be repeated in two weeks.

Allergy Testing

Allergy testing is available including the most commonly found allergens either as a full panel or by using a very cost-effective "screen and expand" approach.

Please see www.allervet.co.uk or contact the laboratory for details.



Oncology





Oncology

Thymidine Kinase (Canine and Feline)

Serum thymidine kinase is a tumour marker for canine and feline malignant lymphoma. Thymidine kinase (TK1) is a cytosolic enzyme associated with cellular proliferation. It markedly increases in activity during the G1 and S phases of the cell cycle. Thymidine kinase enters the circulation from cycling cells or perhaps from cells dying during replication. In essence, then, TK is a marker of the volume of replicating cells in the animal and consequently an indirect measure of neoplastic burden.

The relationship between serum TK and a significant increase in cell replication allows the test to be considered as a diagnostic adjunct for neoplasia such as malignant lymphoma. However, suggestions for its use focus more on prognosis and therapeutic monitoring.

Those lymphoma patients with the highest TK activities are likely to have the highest neoplastic burden or replication rate and these may have a worse prognosis than those lymphoma patients with lower levels.

Thymidine kinase can also be used during the initial phase of chemotherapy to document a decrease from pre-treatment values, indicating the "success" of early therapy and informing the remaining protocol. Frequent monitoring of TK during or between chemotherapeutic cycles should also identify relapse (perhaps pre-clinically) or help differentiate the origin of a clinical deterioration.

Thymidine kinase may not have perfect diagnostic specificity when used as a diagnostic test for lymphoma because of the possible association with other conditions of increased cell turnover including inflammation. However, in the absence of leukocytosis the likelihood of a neoplastic origin is much increased.

Reference:

VonEuler, Einarsson, Olsson, Lagerstedt and Erikson (2004) Serum Thymidine Kinase Activity in Dogs with Malignant Lymphoma: A potent marker for prognosis and monitoring the disease. *Journal of Veterinary Internal Medicine* **18:**696-702

Thymidine Kinase (Equines)

Lymphoma is the most common malignancy in horses and mainly presents as alimentary, cutaneous, mediastinal or multicentric forms.

Thymidine Kinase (TK) has been shown to be significantly higher in horses with lymphoma than clinical normal horses or horses with inflammatory disease or non-haematopoietic neoplasia. A TK level above the reference range appears to be suggestive of the presence of lymphoma in horses but it cannot be ruled out in horses with normal TK values.

Reference:

Larsdotter, Nostell and Von Euler (2015) Serum thymidine kinase activity in clinically healthy and diseased horses: A potential marker for lymphoma. *The Veterinary Journal*



References

A list of references used in the preparation of this manual is available on request.





Notes	
	NationWide LABORATORIES
	01223 493400 www.thehormonelab.co.uk
	www.aienoimonelab.co.ak
MARKET BY	AND



Appendices





Appendices

i. Conversion Factors

Analyte	Per Litre Units	Conversion			Traditional Units
ACTH	pmol/L	divide /L result by	0.2222	to get Trad. units	pg/mL
Aldosterone	pmol/L	н	2.774	11	pg/mL
Calcium	mmol/L	н	0.25	11	mg/dL
Cortisol	nmol/L	11	2.759	11	ng/mL
Cortisol	nmol/L	11	27.59	11	ug/dL
Creatinine	mmol/L	11	0.0884	11	mg/dL
Glucose	mmol/L	11	0.055	11	mg/dL
IGF-1	nmol/L	11	0.131	11	ng/mL
Insulin	pmol/L	11	7.175	11	uIU/mL
Parathyroid hormone	pmol/L	11	0.105	11	pg/mL
Phenobarbital	umol/L	н	4.306	11	ug/mL
Progesterone	nmol/L	н	3.18	11	ng/mL
Testosterone	nmol/L	н	3.467	11	ng/mL
Total T4	nmol/L	11	1.287	11	ng/mL
Total T4	nmol/L	11	12.87	11	ug/dL
Total T3	nmol/L	11	1.536	П	ng/mL
Total T3	nmol/L	11	0.01536	11	ug/dL
Free T4	pmol/L	11	1.287	11	pg/mL

For example, 220 nmol/L Cortisol = $(220 \div 27.59) = 8.0 \text{ ug/dL}$

To convert from traditional to per litre units, multiply traditional result by factor listed above e.g. 8.0 ug/dL Cortisol = $8.0 \times 27.59 = 220 \text{ nmol/L}$

ii. Assay Sample Requirements

Assay	Sample Type	Special Requirements	Volume (ml)
Thyroid			
T4	SH	None	0.3
T3	SH	None	0.5
TSH	SH	None	0.5
FT4ED	SH	None	0.5
Thyroid Antibodies	SH	None	0.2
TGAA	SH	None	0.2

Adrenocortical			
Cortisol (includes ACTH Stim & DEX Screening)	SH	None	0.3
Urinary Cortisol	Urine	None	0.3
ACTH Endogenous (Canine/Feline)	E	Frozen Separated EDTA Plasma	1.0
ACTH Endogenous (Equine)	Е	Separated EDTA Plasma	1.0
OHP	SH	None	0.3
Aldosterone	SH	None	1.0
Renin	E	Frozen Separated EDTA Plasma	2.0

SH	No Gel Tubes	0.5
SH	None	0.6
SH	None	1.0
SH	None	0.6
Serum only	None	0.5
SH	None	0.5
SH	On Ice Packs	0.5
	S H S H S H Serum only S H	S H None S H None S H None Serum only None S H None

Diabetes			
Insulin	S	Starved Sample	0.3
Insulin/Glucose Ratio	S & Fluoride Oxalate*	Starved Sample	S = 0.3 and FO = 0.3
Insulin Antibodies	S	None	0.5
Fructosamine	SH	None	0.3

Gastrointestinal			
Canine TLI	S	Starved Sample	0.4
Feline TLI	S	Starved Sample	0.2
Vitamin B12/Folate	S	None	1.0
Bile Acids	SH	Post sample 2 hrs after fatty meal	0.3
Canine and Feline PLI	S	Starved Sample	0.5

Assay	Sample Type	Special Requirements	Min Sample Volume (ml)
Calcium Disorders			
Parathyroid Hormone (PTH)	E	Frozen EDTA Plasma	1.0
Parathyroid Hormone Related Peptide	E	Frozen EDTA Plasma	1.0
Ionised Calcium	Separated Serum	Anaerobic Collection. Plain tube filled to brim - no air	0.5
25 OH Vitamin D	S	Starved Sample, Separated Serum Stored Frozen	0.5
25 OH Vitamin D and 1.25 OH Vitamin D (Calcitrol)	S	Starved Sample, Separated Serum Stored Frozen"	2.0

Serology / Immunology			
Acetylcholine Receptor Antibody	S	None	0.5
2M Antibodies	S	None	0.5
Sarcoptic Mange Antibody	S	None	0.2

Therapeutic Drug Monitoring			
Phenobarbitone	S	No Gel Tubes	0.3
Potassium Bromide (KBr)	SH	No Gel Tubes	0.3
Levetiracetam	S	No Gel Tubes	0.5
Digoxin	S	No Gel Tubes	0.5
Cyclosporin	E (Whole Blood)	Do not post Friday	1.0

Miscellaneous			
IGF-1	SH	None	0.2
Erythropoietin	S	Frozen Serum	1.0
Thymidine Kinase	S	None	0.3
Gastrin	S	Frozen Separated Serum	1.0
Normetanephrine (Urine)	Acidified urine	None	20

Biochemistry Profiles			
All Profiles	SH	None	1.0

Sample Guide

E = EDTA
S = Serum
H = Heparin

An ideal sample volume is 1ml of separated serum or plasma.

This allows for assays to be repeated if necessary or further analysis to be added if required.

*A serum sample can also be submitted for Glucose if this has been separated within one hour of collection



iii. Specialist Biochemistry Profiles

NationWide Specialist Laboratories



Biochemistry Profiles

Canine Addisonian ACTH Stimulation Test, Aldosterone Pre and Post, Sodium, Potassium, Na:K Ratio **Profile - ADDP** Senior Canine Cushing's ACTH Stimulation Test, ALT, ALP, Total Protein, Albumin, Globulin, **Profile - CUSH** A:G Ratio, Urea, Creatinine, Glucose, Calcium, Phosphorus, Cholesterol Total T4, ALT, ALP, Total Protein, Albumin, Globulin, A:G Ratio, Urea, Senior Canine Profile - BK9P Creatinine, Glucose, Calcium, Phosphorus, Cholesterol Canine/Feline Pancreatic Profile Basic - PANC Amylase, Lipase, Calcium, Albumin Feline ACTH Stimulation Test, Aldosterone on basal Feline Addisonian Profile - BFAP sample, Sodium, Potassium, Na:K Ratio Total T4, ALT, ALP, Total Protein, Urea, Feline Hyperthyroid Profile - Silver - BFSP Creatinine, Phosphorus Total T4, Free T4ED, ALT, ALP, Total Protein, Feline Hyperthyroid Profile - Gold - BFGP Urea, Creatinine, Phosphorus Feline Pancreatic Profile Plus - PANP Feline TLI, Amylase, Lipase, Calcium, Albumin **Senior Feline** Total T4, ALT, ALP, Total Protein, Albumin, Globulin, A:G Ratio, Urea, Creatinine, Glucose Calcium, Phosphorus, Sodium, Potassium, Na:K Ratio, Fructosamine, Billirubin **Profile - BFSR** Pre Op Profile - BPOP ALT, Urea, Creatinine, Bile Acid **Small Animal Kidney** Total Protein, Albumin, Globulin, A:G Ratio, Urea, Creatinine, Glucose, Calcium, Phosphorus, Sodium, Potassium, Na:K Ratio. **Profile - BKP Small Animal Biochemis-**ALT, ALP, Total Protein, Albumin, Globulin, A:G Ratio, Urea, Creatinine, Glucose, Calcium, Phosphorus, Cholesterol, Bilirubin try Profile - BSAP Urea, Creatinine, Phosphorus, Total Protein, Albumin, Globulin, **Equine Kidney Profile - EKP** Calcium, Sodium, Potassium, Chloride AP, AST, GGT, Bilirubin, GLDH, LDH, Bile Acid, Total Protein, **Equine Liver Profile - ELP** Albumin, Globulin, Fibrinogen, Serum Amyloid A

AP, AST, GGT, Bilirubin, GLDH, Total Protein, Albumin, Globulin, Equine Basic Profile - EGP Urea, Creatinine, Phosphorus, Calcium, CK, Bile Acid, LDH.

Equine Basic Profile AP, AST, GGT, Bilirubin, GLDH, Total Protein, Albumin, Globulin, Urea, Creatinine, Phosphorus, Calcium, CK, Bile Acid, LDH, Fibrinogen, Serum Amyloid A. Plus - ESP

Equine Laminitis Profile Plus - LAMB

ACTH, Insulin, Cortisol, Glucose, Urea, Creatinine, Calcium, Phosphorus, Sodium, Potassium, Total Protein, Albumin

For further information contact **NationWide Specialist Laboratories** Unit 2, Sawston Park, London Road, Pampisford, Cambridge, CB22 3EE T. 01223 493400 email: csls.info@nwlabs.co.uk or visit www.thehormonelab.com

For all these profiles, the Total T4, Free T4 by equilibrium dialysis, Cortisol, Aldosterone and ACTH assays will be run using RIA (radioimmunoassay) methods. We believe this to be the best method to accurately measure hormones.

iv. Laboratory Request Form

	Owners name: Animal Name / Reference: Species: Breed:			NationWide
	Age: Sex:	☐ Neut		LABORATORIES
M –	Clinical Details / Current Therapy:		NationWide Spe Assay Request	ecialist Laboratories Form
	Sample: Serum Plasma Date: Time(s):	☐ Urine	Clinician: Practice: Address:	
	SUPPLIES REQUESTED:		E-mail for results	s:
	AYROID CANINE PROFILES Platinum Profile (T4, TSH, FT4ED, TGAA) Copper Profile (FT4ED, TSH, T4) Gold Profile (TSH, FT4ED, TGAA) Silver Profile (TSH, FT4ED) Bronze Profile (TSH, FT4ED) Bronze Profile (TGAA, FT4ED) Bronze Profile (TGAA, FT4ED) Greyhound Thyroid Investigation (T4, T3, TSH, TGAA) Greyhound Thyroid Monitoring (T4, T3, TSH, TGAA) Greyhound Thyroid Profile - Silver (Total T4, ALT, ALP, Total Protein, Urea, Creatinine, Phosphorus) Feline Hyperthyroid Profile - Gold (As Silver plus FT4ED) Feline Hyperthyroid Profile Basic (T4, FT4ED) INDIVIDUAL TESTS Thyroxine (T4) Free T4 by dialysis (FT4ED) TSH TRH/TSH stimulation test TGAA Anti T4 Antibodies Anti T3 Antibodies THAB (Anti T4 Antibodies & Anti T3 Antibodies) Triiodothyronine (T3) DRENOCORTICAL Canine ACTH stim Test (0, 1 & 3 hrs) Dexamethasone screening test (0, 3 & 8 hrs) ACTH (EDTA only, special transport) Feline combined HDD/ACTH test Cortisol/Creatinine Ratio (urine) Cortisol Saliva Cortisol SHAP profile (Cort x 2, OHP x 2) Aldosterone Aldosterone stim test (pre & post) Hypertension/Hypokalaemia Investigation (Plasma Renin, Aldosterone, Na:K) Renin (EDTA plasma/special transport) Canine Addisonian Profile (2 x Cort/Aldo, Na:K) Ferret Adrenal Profile (E2, OHP)	Testosterone hCG stimulation te AMH - Canine/Feli Relaxin - Canine o THERAPEUTIC MONI (Do not use serum gel Phenobarbitone Phenobarbitone Phenobarbitone, GLDH, Bile Acid) Potassium Bromid Vitamin Bromid Urea, Creatinine, F Digoxin Cyclosporin Levetiracetam GASTROINTESTINAL Canine TLI (serum) Vitamin B12/Folate Feline TLI (serum) Canine/Feline PLI Bile Acids (pre & p Gastrin (special sa Lipase DIABETES Insulin Insulin/Glucose Ra Insulin Antibodies Fructosamine CALCIUM HOMEOST PTH PTH/PTHrP Ionised Calcium (IC PTH/ICa PTH/PTH/PICa 25-OH-Vitamin D 25-OH-Vitamin D SEROLOGY/IMMUNO ACR antibodies Sarcoptic Mange a Canine Prostate Sp	test (2 x E2) Feline hcG Stim test) set (2 x Testosterone) ne nily TORING tubes) conitoring Profile ALP, ALT, Albumin, GGT, e e Monitoring Profile (KBR, Phosphorous) c (serum) (serum) ost) mpling requirements) attio ASIS (Special Sampling) classical sampling of the service of the ser	EQUINE SPECIALITIES Endogenous ACTH – Equine (EDTA plasma) Laminitis Profile (ACTH, Insulin, Cortisol & Glucose) Laminitis Profile Plus (as above + chemistries) Equine O/N DXM Test (2 x Cort 0 & 22hrs) Oestrone Sulphate Rig hCG Stimulation (2 x Testosterone) Oestrone Sulphate Pregnancy PMSG Inhibin AMH (Anti-Müllerian Hormone) GCTA Profile (Progesterone, AMH) GCTB Profile (AMH, Progesterone & Testosterone) GCT Profile Plus (Progesterone, Testosterone) GCT Profile Plus (Progesterone, Testosterone, Inhibin & AMH) MISCELLANEOUS IGF1 (serum)

Units S5 and S6, The Saxon Building, Eastern Counties Leather Industrial Estate, London Road, Pampisford, Cambridge, CB22 3FJ.

Tel: 01223 493400 • Email: csls.info@nwlabs.co.uk

NationWide Laboratories is a trading business of National Veterinary Services Ltd. a reg. Co. in England & Wales (Reg No. 8489997). Reg. Add: Unit 4, Jamage Ind. Est., Talke Pits, Stoke-on-Trent, Staffs, ST7 1XW. VAT No.GB 169 1627 83

v. Laboratory Assay Schedule

Monday	Tuesday	Wednesday	Thursday	Friday
T4	T4	T4	T4	T4
Cortisol	Cortisol	Cortisol	Cortisol	Cortisol
Insulin	TSH	OHP	Insulin	TSH
Canine ACTH	Feline TLI	Sarcoptic Mange	Feline TLI	Т3
		PMSG	CPSE	T3/T4 Antibodies
All Biochemistries	All Biochemistries	All Biochemistries	All Biochemistries	All Biochemistries
Ionised Calcium	Ionised Calcium	Ionised Calcium	Ionised Calcium	Ionised Calcium
		AMH		Relaxin
				Progesterone*
Assays Set up Today Include	Assays Set up Today Include	Assays Set up Today Include	Assays Set up Today Include	Assays Set up Today Include
Oestrone Sulphate	Testosterone	Equine ACTH	TGAA	IGF1 (to incubate over weekend)
FT4ED	Oestradiol		FT4ED	
	Aldosterone	PTHrP		
PTH	РТН	PTH	PTH	РТН
(Androstenedione run on alternate Mondays) - please check with Customer Services for details				
Results Available PM	Results Available PM	Results Available PM	Results Available PM	Results Available PM
IGF1	Oestrone Sulphate	Testosterone	Equine ACTH	TGAA
	FT4ED	Oestradiol	PTHrP	FT4ED
		Aldosterone		

Please Note

Should a result be needed urgently please contact Customer Services as it may be possible to run the assay sooner than the above schedule

This service will incur an extra charge (usually double list price)

*Premate Progesterone Assayed on day of receipt

Notes	
	NationWide LABORATORIES
	01223 493400 www.thehormonelab.co.uk
CONTRACTOR OF THE SECOND SECON	
AND STALL	NA JOR

Notes	®
	NationWide LABORATORIES
	01223 493400
	www.thehormonelab.co.uk
-	
4	
	100
	9/11/11
KYL by	A DAY OF THE

SPECIALIST LABORATORY **SERVICES INFORMATION MANUAL**

20th EDITION



Telephone:

+44 (0)1223 493400

Fmail:

csls.info@nwlabs.co.uk

Website:

www.thehormonelab.com

Address

NationWide Specialist Laboratories Units S5 & S6, The Saxon Building, Eastern Counties Leather Ind Est. London Rd, Pampisford, Cambs, CB22 3FJ.

Opening times

Monday to Friday: 8am - 4.30pm

Additional services from NWL



NationWide Laboratories (Head Office)

Tel: 01253 899215 Fax: 01253 891934 Lancefield House, 23 Mains Lane Little Singleton, Poulton-le-Fylde Lancashire FY6 7LJ

info@nwlabs.co.uk Fmail:

Saturday:

Monday to Friday: 8.30am to 5.30pm 8.00am to 12.00pm

nwlabs.co.uk WWW.

A telephone answering service is in operation when the laboratories are closed.